


FFAP Test Mix

Column: Zebron™ ZB-FFAP, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea
Phase: Nitroterephthalic Acid Modified Polyethylene Glycol
Dimensions: 30 meters x 0.32 mm x 0.25 µm
Order No: 7HM-G009-11
Oven Profile:

Carrier Gas: Hydrogen, 1.93 mL/min
Injection: Split 50.0000000000:1 0 µL @ 250°C
Detection: Refractive Index (255°C)



Products used in this application:

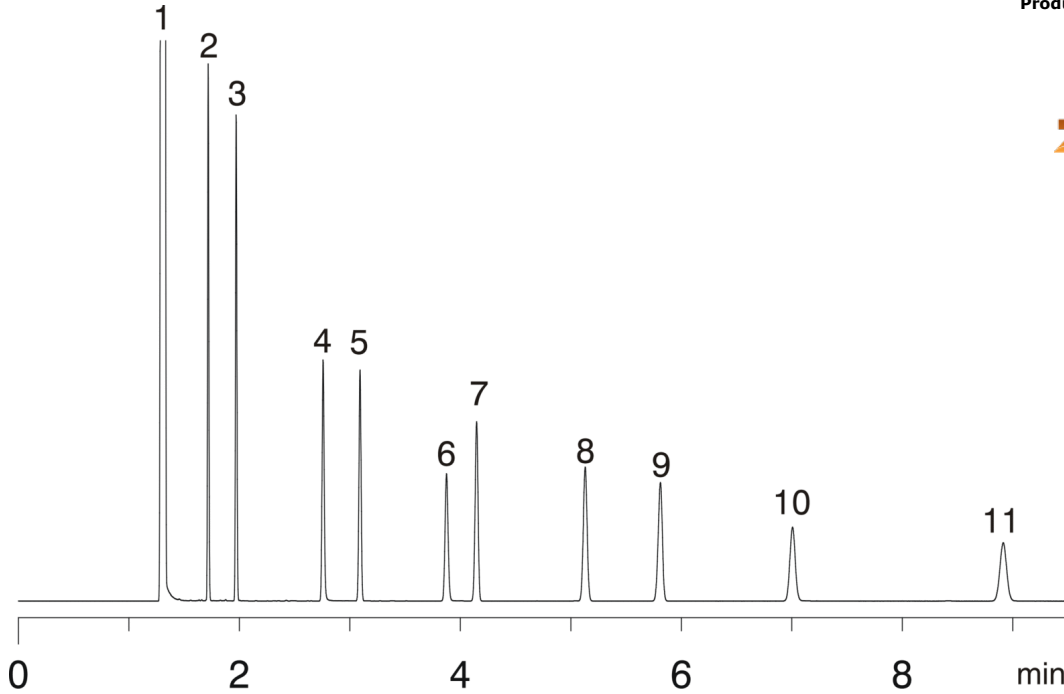

FFAP Test Mix

ANALYTES:

- 1 Methane
- 2 Tridecane
- 3 Tetradecane
- 4 1-Octanol
- 5 Methyl decanoate
- 6 Isovaleric acid
- 7 Methyl undecanoate
- 8 Naphthalene
- 9 Methyl dodecanoate
- 10 1-Hexanoic acid
- 11 2,6-Dimethylphenol

