

## Increase Sample Through-Put (fig. 2)

**Column:** Luna® 3 µm C18(2) 100 Å, LC Column 20 x 4 mm, Ea

**Dimensions:** 20 x 4 mm ID

**Order No:** 00M-4251-D0

**Elution Type:** Isocratic

**Eluent A:** Water with 0.1% ammonium hydroxide, pH 10

**Eluent B:** Acetonitrile with 0.1% ammonium hydroxide, pH 10

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 0        | 80         | 20    |       |

**Flow Rate:** 2 mL/min

**Col. Temp.:** 25 °C

**Detection:** Refractive Index @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

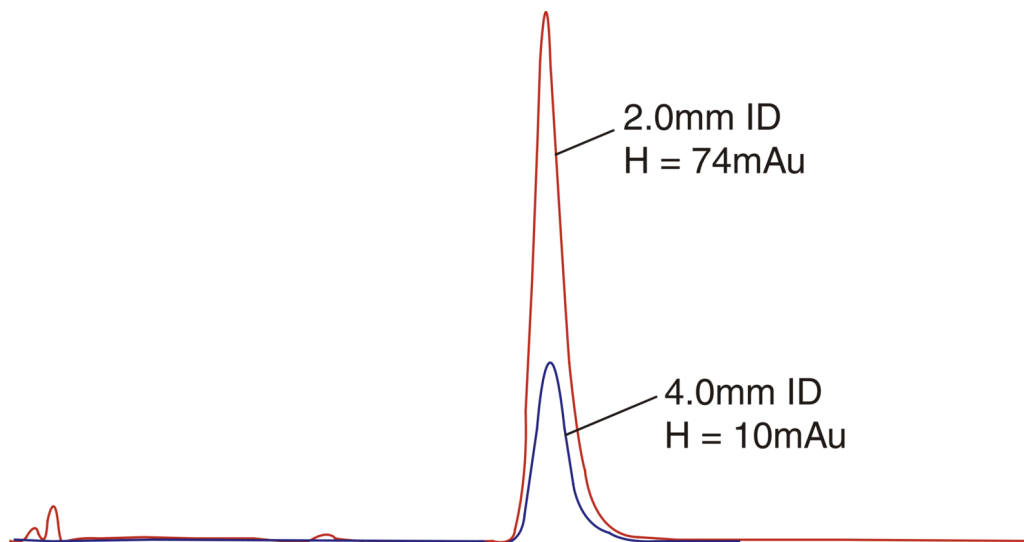

### ANALYTES:

1 Pindolol

