

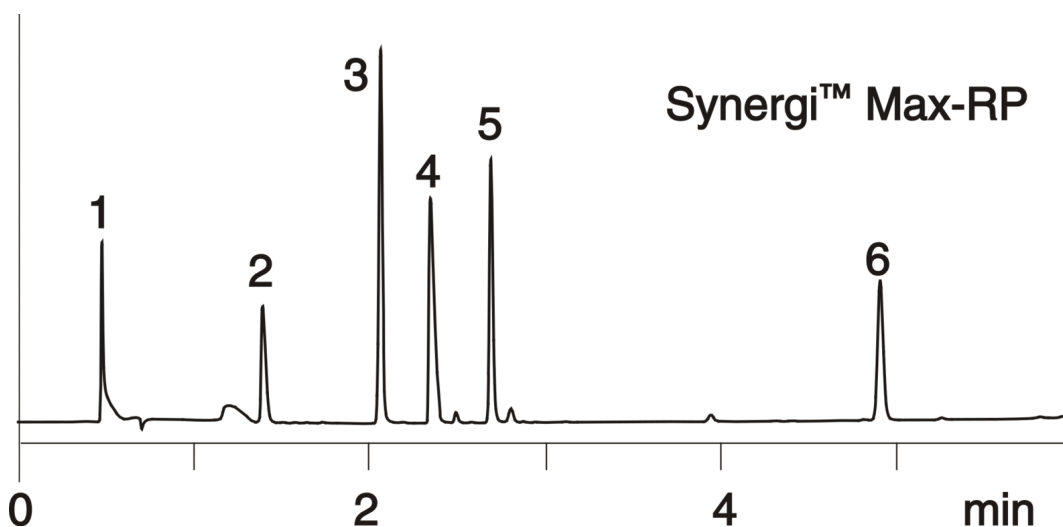
Acidic and basic compounds using gradient**Column:** SynergiTM 4 μ m Max-RP 80 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4337-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% Formic acid**Eluent B:** Acetonitrile with 0.1% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	0	100

Flow Rate: 1.5 mL/min**Col. Temp.:** 30 °C**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Thiourea
- 2 Codeine
- 3 Chlorpheniramine
- 4 Propranolol
- 5 Desipramine
- 6 Ibuprofen

