


Volatile Organic Compounds**Column:** Zebron[™] ZB-624, GC Cap. Column 30 m x 0.25 mm x 1.40 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 1.4 µm**Order No:** 7HG-G005-27**Oven Profile:** 32Å°C for 4 min 80Å°C for 12 min 100Å°C for 22 min**Carrier Gas:** Constant Flow Helium, 0.88 mL/min**Injection:** Split 200.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (300°C)**Analyst Note:** Sample analytes 1000ppm EPA Method 624 mixture**Products used in this application:**

Volatile Organic Compounds

ANALYTES:

- 1 1,1-Dichloroethene
- 2 Methylene chloride
- 3 1,1-Dichloroethane
- 4 Carbon tetrachloride
- 5 Benzene
- 6 Trichloroethene
- 7 1,2-Dichloropropane
- 8 2-Chloroethyl vinyl ether
- 9 1,1,2-Trichloroethane
- 10 Tetrachloroethane
- 11 Dibromochloromethane
- 12 Chlorobenzene

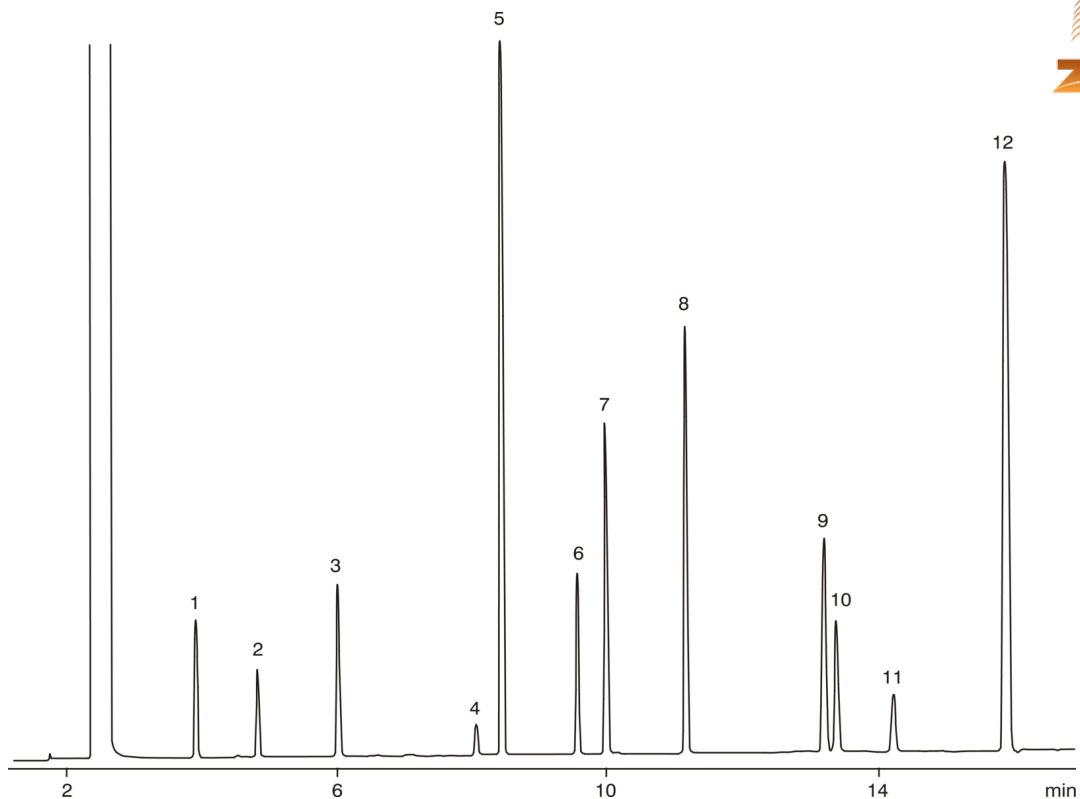