


Pesticides: EPA Method 8081A**Column:** Zebron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G004-11**Oven Profile:** 270Â°C for 38 min 270Â°C for 40 min**Carrier Gas:** Constant Flow Helium, 1.09 mL/min**Injection:** Split 20.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (25°C)**Products used in this application:**

Pesticides: EPA Method 8081A

ANALYTES:

- 1 alpha-BHC
- 2 gamma-BHC
- 3 beta-BHC
- 4 Heptachlor
- 5 delta-BHC
- 6 Aldrin
- 7 Heptachlor epoxide
- 8 Endosulfan I
- 9 4,4'-DDE
- 10 Dieldrin
- 11 Endrin
- 12 4,4'-DDD
- 13 Endosulfan II
- 14 4,4'-DDT
- 15 Endrin aldehyde
- 16 Endosulfan sulfate
- 17 Methoxychlor
- 18 Endrin ketone

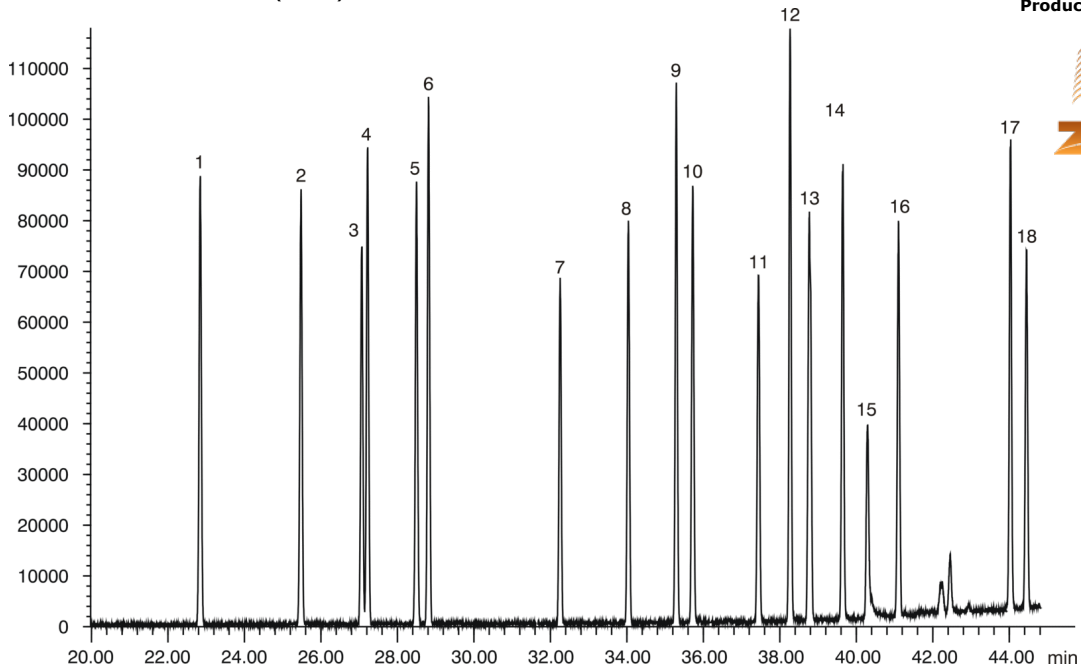