

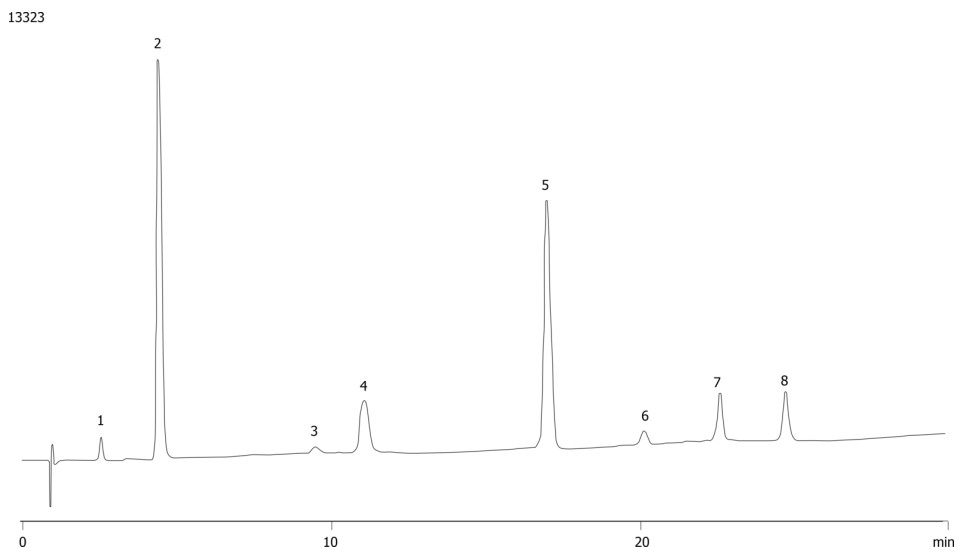
Phenolic compounds**Column:** Prodigy[™] 5 µm ODS-2 150 Å, LC Column 75 x 4.6 mm, Ea**Dimensions:** 75 x 4.6 mm ID**Order No:** 00C-3300-E0**Elution Type:** Gradient**Eluent A:** Acetonitrile**Eluent B:** Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	25	75
	2	11	30	70
	3	30	60	40
	4	30.1	100	0
	5	45	100	0

Flow Rate: 1 mL/min**Col. Temp.:** 50 °C**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**Products used in this application:****Prodigy[™]****ANALYTES:**

- 1 Phenol
- 2 o-Cresol
- 3 IPP
- 4 p,p-BPA
- 5 o,p-BPA
- 6 IPP2
- 7 BPX
- 8 COD

