

## Brompheniramine

**Column:** Chirex® (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3022-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol-TFA 20:1  
**Gradient**  


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

**CHIREX**®

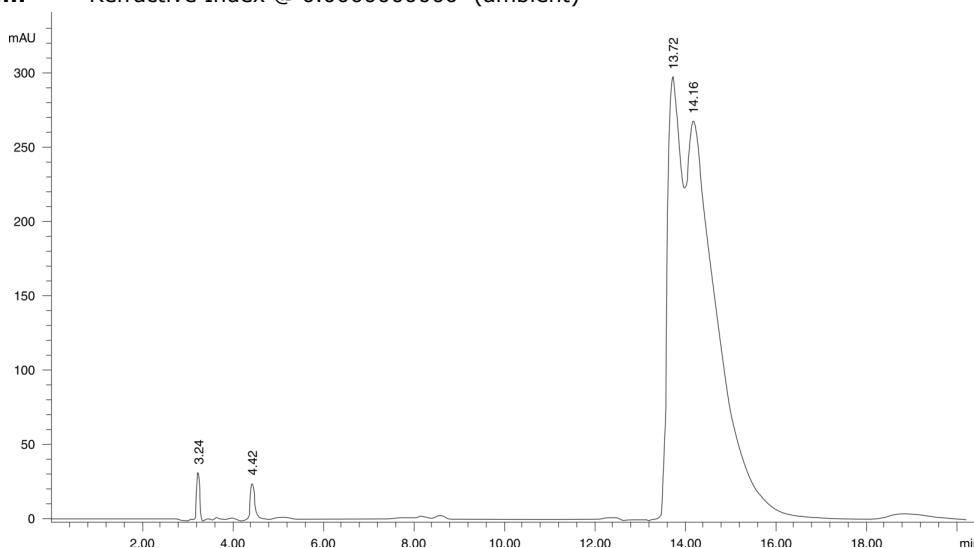

### ANALYTES:

- 1 Brompheniramine  
Retention Time: 13.72 min
- 2 Brompheniramine  
Retention Time: 14.16 min

