

## Diperodon

**Column:** Chirex<sup>®</sup> (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3014-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA 60:35:5 ethanol-TFA 20:1  
**Gradient**


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 0.7 mL/min  
**Col. Temp.:** ambient  
**Detection:** Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

**CHIREX<sup>®</sup>**

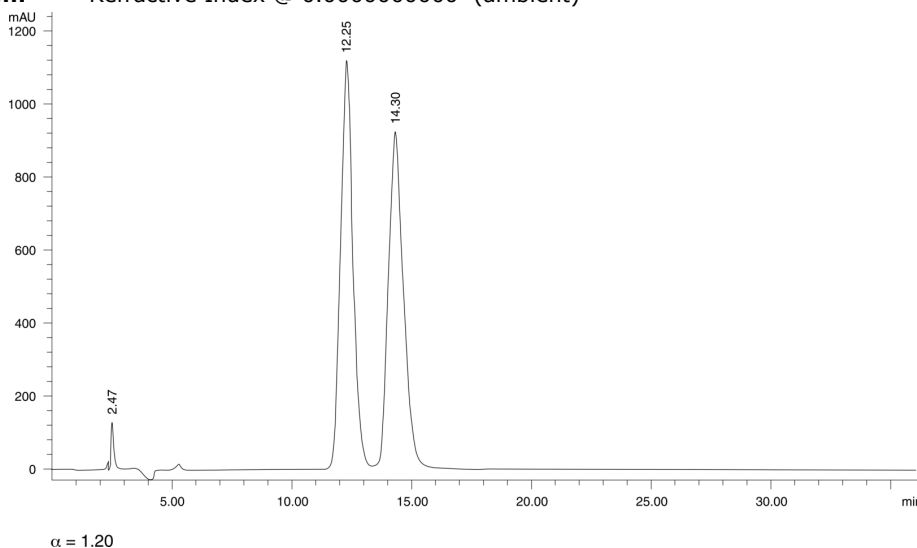

### ANALYTES:

- 1** Diperodon  
Retention Time: 12.25 min
- 2** Diperodon  
Retention Time: 14.3 min

