


Diperodon


Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3020-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA 77:20:3 ethanol-TFA 20:1
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)

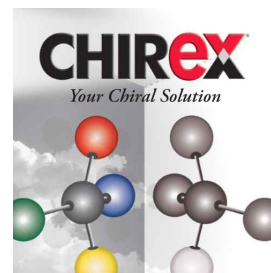

Products used in this application:

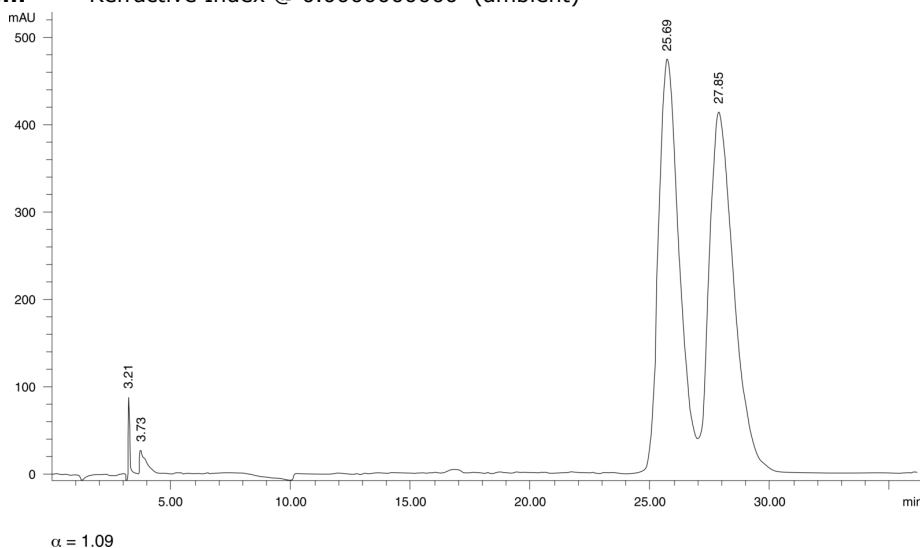CHIREX®

ANALYTES:

- 1 Diperodon
Retention Time: 25.69 min
- 2 Diperodon
Retention Time: 27.85 min

