


## Metaproterenol

**Column:** Chirex<sup>®</sup> (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3022-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol-TFA 20:1  
**Gradient**


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Refractive Index @ 0.0000000000 (ambient)

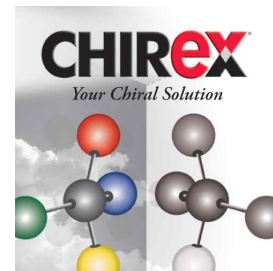

Products used in this application:

**CHIREX<sup>®</sup>**

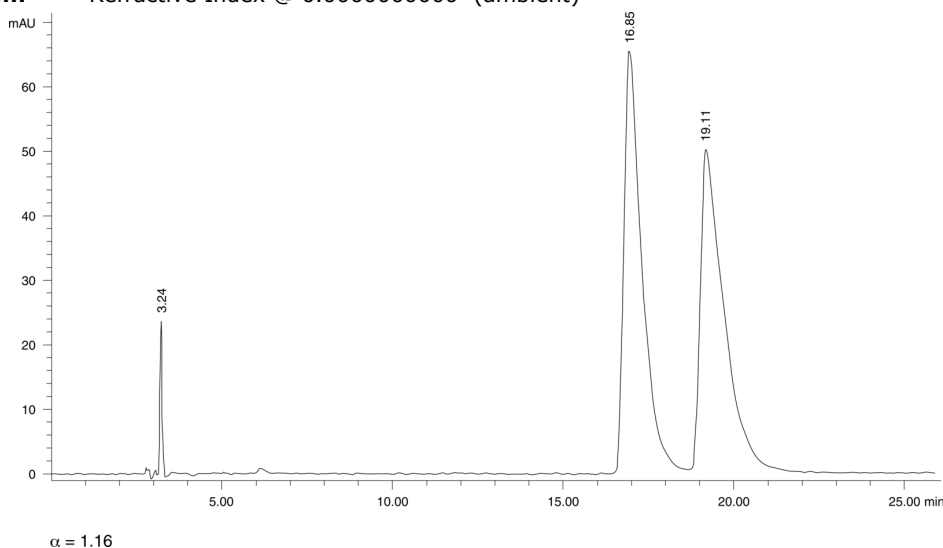

### ANALYTES:

- 1 Metaproterenol  
Retention Time: 3.24 min
- 2 Metaproterenol  
Retention Time: 16.85 min
- 3 Metaproterenol  
Retention Time: 19.11 min

