

HPLC Application


ID No.: 13662

Metoprolol

Column: Chirex® (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3022-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1
Gradient


Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)

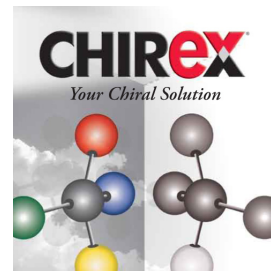

Products used in this application:

CHIREX

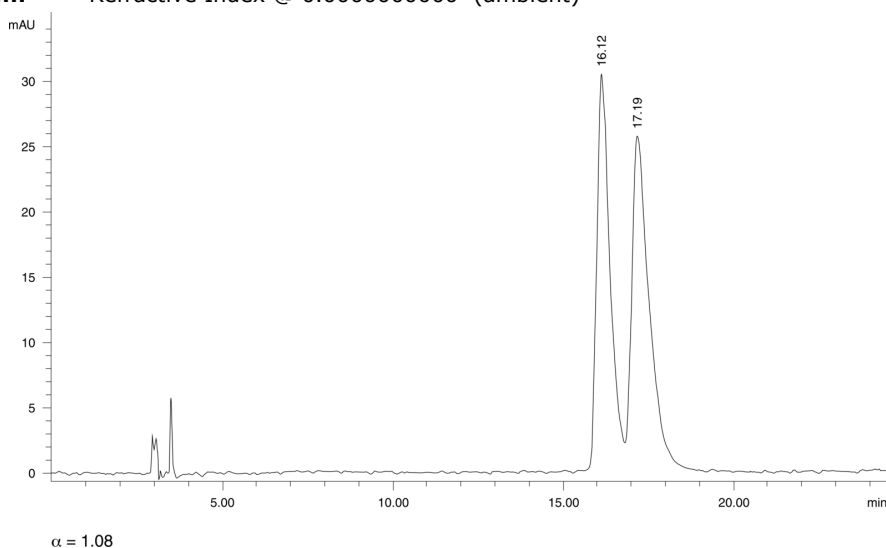

ANALYTES:

- 1 Metoprolol
Retention Time: 16.12 min
- 2 Metoprolol
Retention Time: 17.19 min

