


Nefopam


Column: Chirex® (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3022-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (58:35:7) w/ethanol-TFA (20:1)
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)

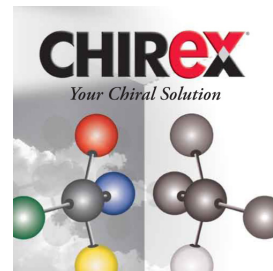

Products used in this application:

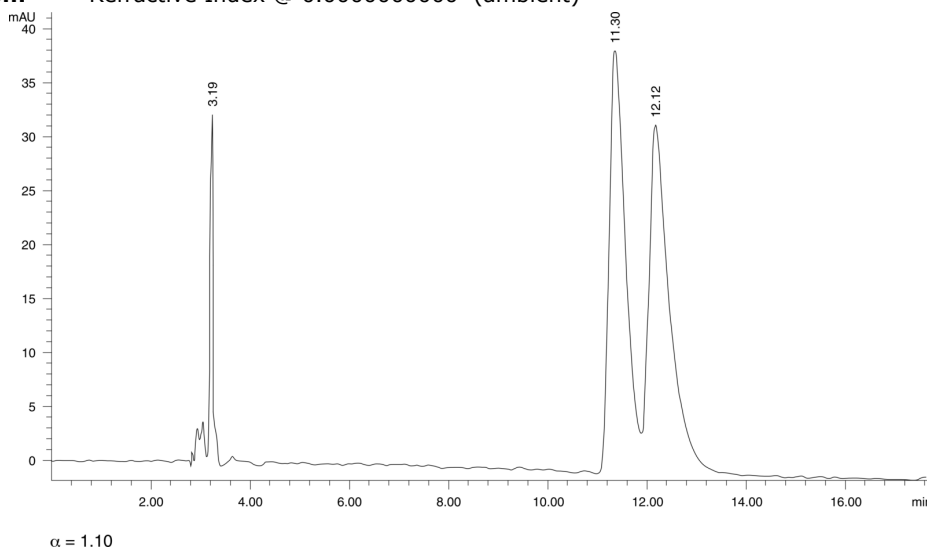CHIREX®

ANALYTES:

- 1 Nefopam
Retention Time: 3.19 min
- 2 Nefopam
Retention Time: 11.3 min
- 3 Nefopam
Retention Time: 12.12 min

