

Tetrahydrozoline

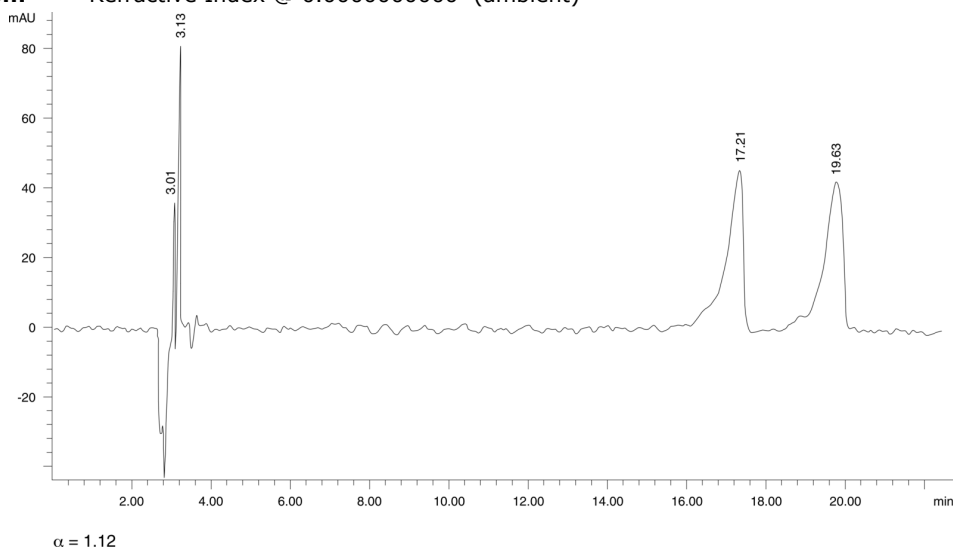
Column: Chirex® (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3022-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (50:35:150 w/ ethanol-TFA (20:1)
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

CHIREX®

ANALYTES:

- 1 Tetrahydrozoline
Retention Time: 17.21 min
- 2 Tetrahydrozoline
Retention Time: 19.63 min

