

HPLC Application

ID No.: 13720

Warfarin

Column: Chirex[®] (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3022-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (88:10:2) ethanol-TFA 20:1
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

CHIREX[®]

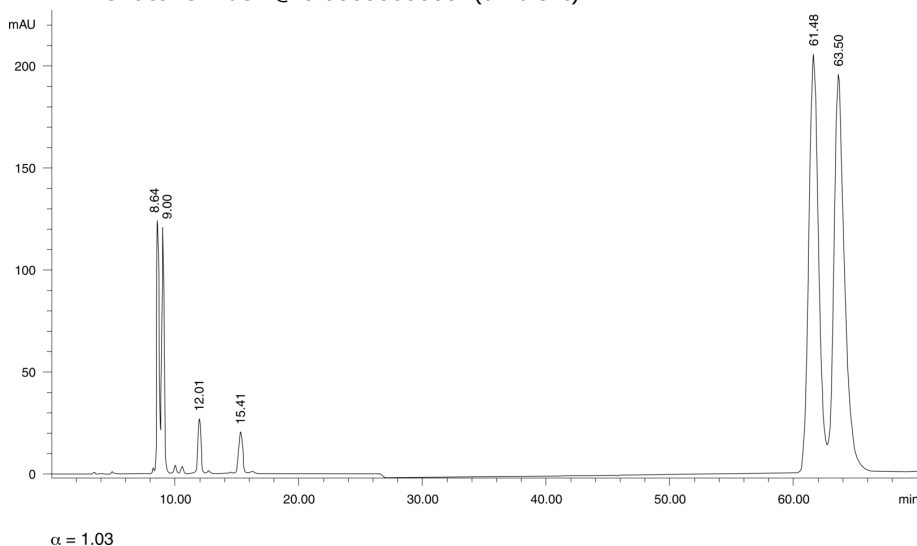

ANALYTES:

- 1 Warfarin
Retention Time: 61.48 min
- 2 Warfarin
Retention Time: 63.5 min

