

Octopamine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 5 mM Copper (II) acetate in water [pH 6.0]



Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

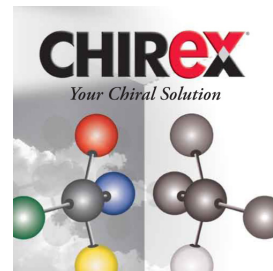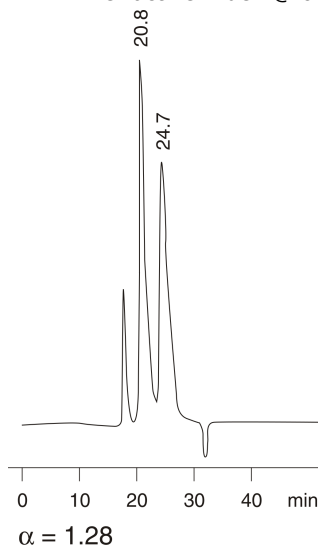

Products used in this application:

CHIREX

ANALYTES:

1 Octopamine

Retention Time: 20.84 min

2 Octopamine

Retention Time: 24.68 min

