

## 2-Amino-n-butyric acid

**Column:** Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-3126-E0

**Elution Type:** Isocratic

**Eluent A:** 1 mM Copper (II) sulfate in water

**Gradient**


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**


**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)

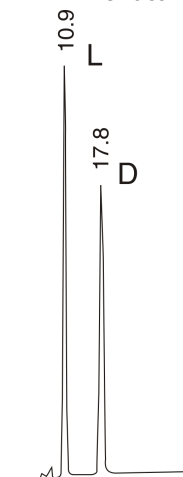

$$\alpha = 1.80$$

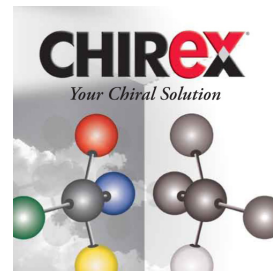

Products used in this application:

**CHIREX**

### ANALYTES:

1 L-2-Amino-n-butyric acid

Retention Time: 10.85 min

2 2-Amino-n-butyric acid

Retention Time: 17.703 min

