

Ornithine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) sulfate in water


Gradient	Step No.	Time (min)	Pct A
----------	----------	------------	-------

Profile:	1	0	100
----------	---	---	-----

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

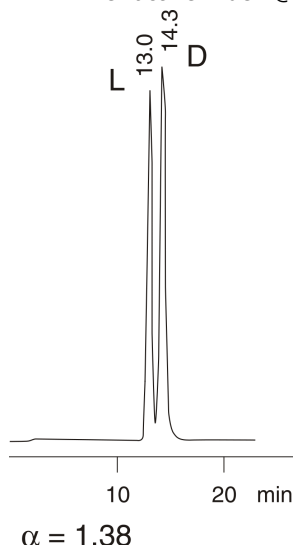

Products used in this application:

CHIREX®

ANALYTES:

1 L-Ornithine

Retention Time: 13.833 min

2 D-Ornithine

Retention Time: 14.283 min

