

N-Acetylvaline

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-3126-E0
Elution Type: Isocratic
Eluent A: 2 mM Copper (II) sulfate in water
Eluent B: Acetonitrile
Gradient


Step No.	Time (min)	Pct A	Pct B
1	0	85	15

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

CHIREX

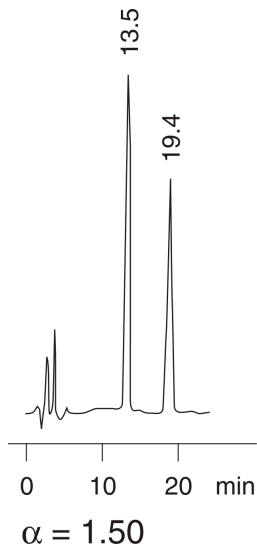

ANALYTES:

- 1 N-Acetylvaline
Retention Time: 13.54 min
- 2 N-Acetylvaline
Retention Time: 19.42 min

