

3-Amino-e-caprolactam

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) sulfate in water

Gradient



Step No.	Time (min)	Pct A
----------	------------	-------

Profile:	1	0	100
-----------------	----------	---	-----

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

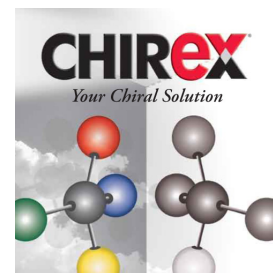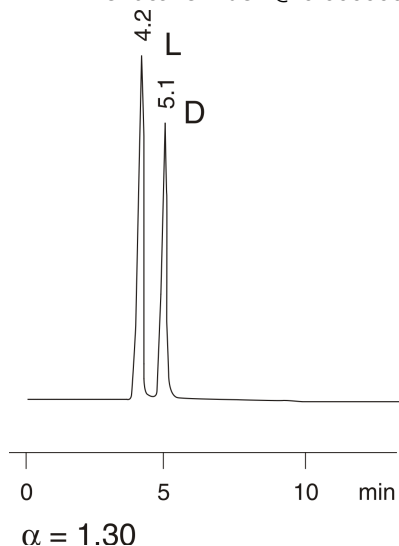

Products used in this application:

CHIREX

ANALYTES:

- 1** L-3-Amino-e-caprolactam

Retention Time: 4.21 min

- 2** Amino-e-caprolactam-d3

Retention Time: 5.05 min

