

Tetrahydro-2-furancarboxylic acid

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) sulfate in water

Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

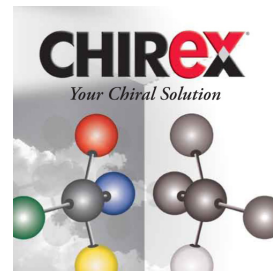

Products used in this application:

CHIREX

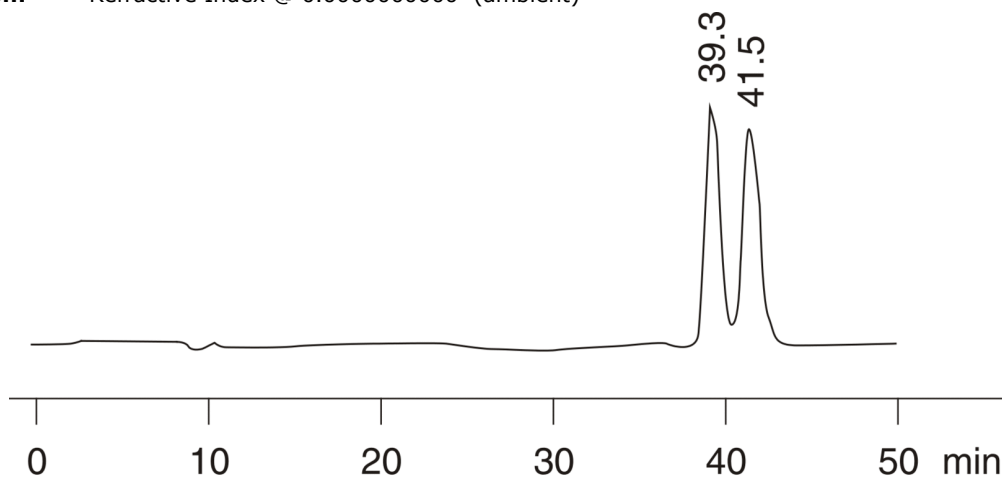

$$\alpha = 1.06$$

ANALYTES:

- 1 Tetrahydro-2-furancarboxylic acid

Retention Time: 39.263 min

- 2 Tetrahydro-2-furancarboxylic acid

Retention Time: 41.483 min

