

HPLC Application

ID No.: 14117

Pantoic acid

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 2 mM Copper (II) sulfate in water


Eluent B: Isopropanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

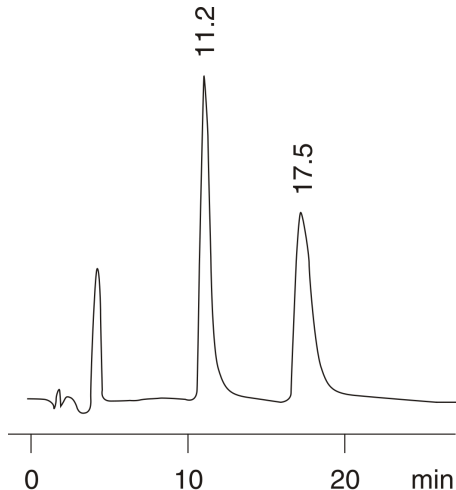

$\alpha = 1.63$

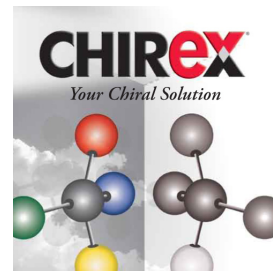

Products used in this application:

CHIREx

ANALYTES:

- 1 Pantoic acid
Retention Time: 11.198 min
- 2 Pantoic acid
Retention Time: 17.485 min