

### Simple Sugars (#3) on Luna Amino

**Column:** Luna® 5 µm NH2 100 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4378-E0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** Water

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 80    | 20    |

**Flow Rate:** 3 mL/min

**Col. Temp.:** 40 °C

**Detection:** Refractive Index @ 0.0000000000 (ambient)

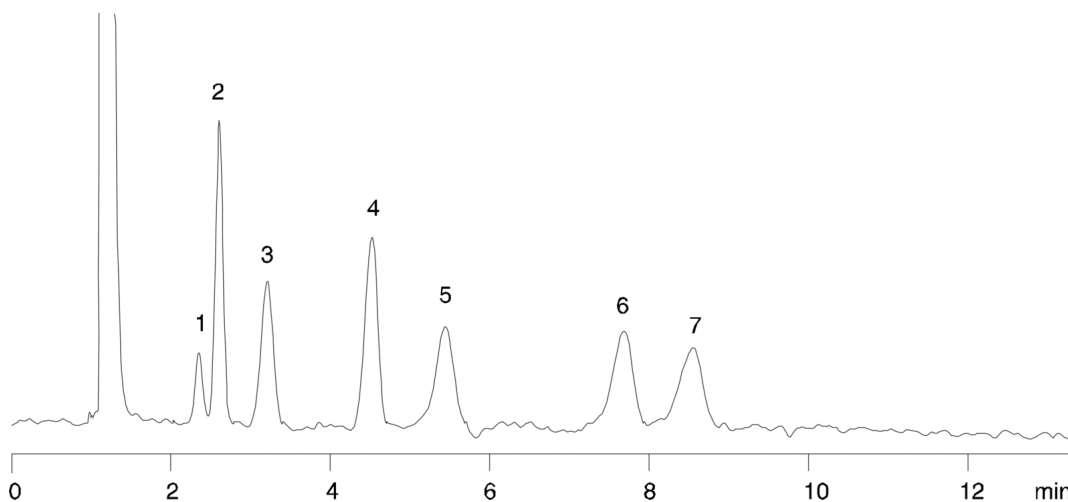
**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, NH2 4 x 3.0mm, 10/Pk Part No.: AJ0-4302

- Holder Part No.: KJ0-4282



**Products used in this application:**

**ANALYTES:**

- 1 Xylose
- 2 Fructose
- 3 Glucose
- 4 Sucrose
- 5 Maltose
- 6 Melezitose
- 7 Raffinose