

Peptide Comparison - Jupiter Proteo

Column: Jupiter® 4 µm Proteo 90 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4396-E0

Elution Type: Gradient

Eluent A: 0.01%TFA in Water

Eluent B: 0.008% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	20	60	40

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282



Products used in this application:

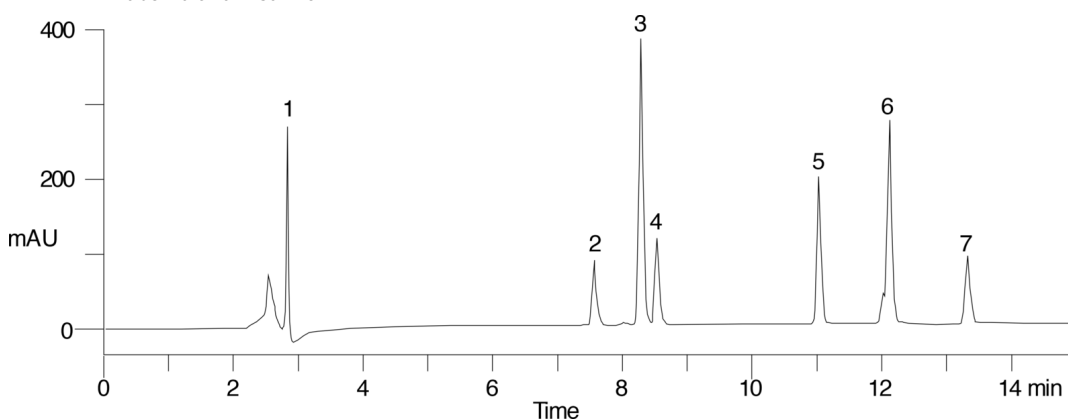

ANALYTES:

- 1 Angiotensinogen fragment 11-14 (Val-Ile-His-Asn)
- 2 [sar1, val5, ala8] Angiotensin II
- 3 [sar1, thr8]-Angiotensin II
- 4 [sar1, ala8]-Angiotensin II
- 5 [asn1, val5]-Angiotensin II
- 6 [sar1]-Angiotensin II
- 7 Angiotensin I

