

TFA Concentration Comparison (3) - Tryptic Digest of Beta-Lactoglobulin


Column: Jupiter® 4 µm Proteo 90 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4396-E0**Elution Type:** Gradient**Eluent A:** 10mM TFA**Eluent B:** 10mM TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	95	5
	3	60	55	45


Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Refractive Index @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

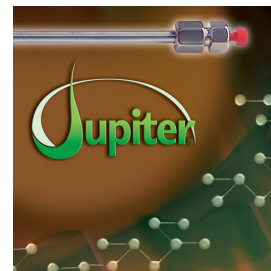

Products used in this application:

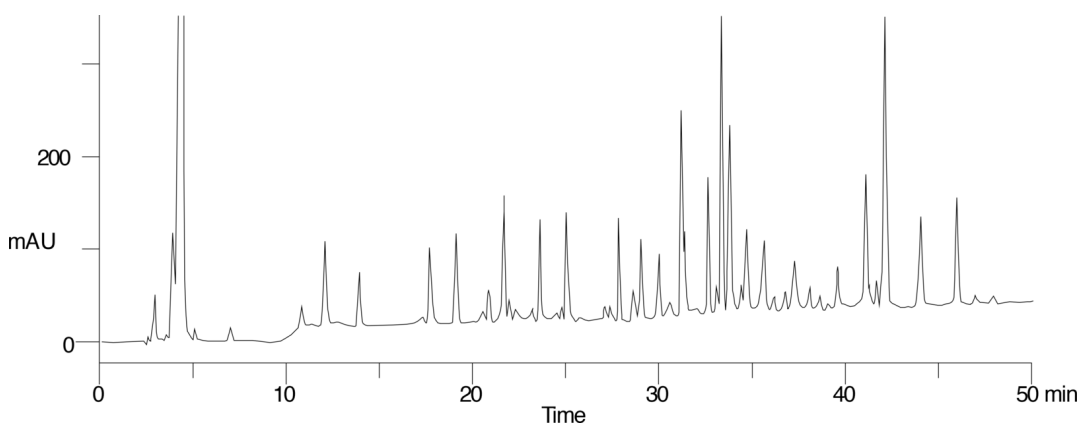

ANALYTES:

1 beta-Lactoglobulin

