

HPLC Application

ID No.: 14478

EPA Method 610, PAH's

Column: SynergiTM 4 μ m Max-RP 80 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4337-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	5	60	40
	3	30	0	100

Flow Rate: 1.55 mL/min


Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

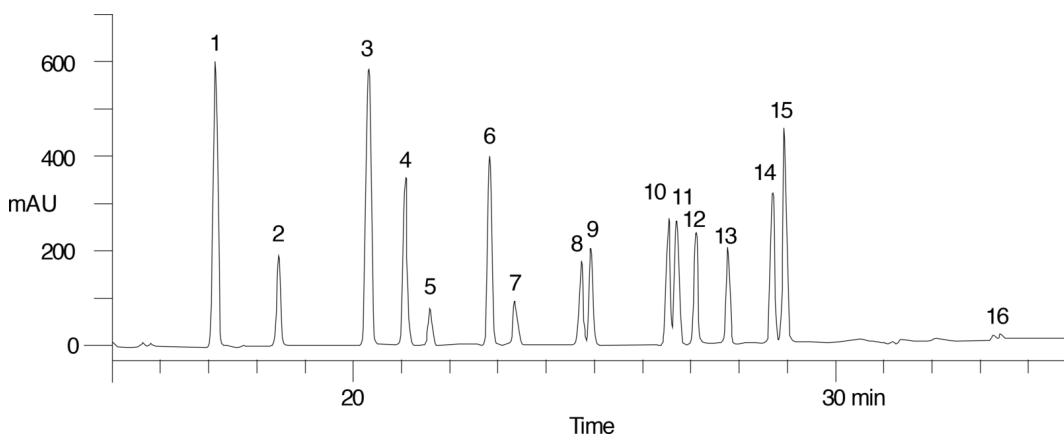

EPA Method 610, PAH's

ANALYTES:

- 1** Naphthalene
- 2** Acenaphthylene
- 3** Acenaphthene
- 4** Fluorene
- 5** Phenanthrene
- 6** Anthracene
- 7** Fluoranthene
- 8** Pyrene
- 9** Benz[a]anthracene
- 10** Chrysene
- 11** Benzo[b]fluoranthene
- 12** Benzo[k]fluoranthene
- 13** Benzo[a]pyrene
- 14** Dibenzo[a,h]anthracene
- 15** Benzo[g,h,i]perylene
- 16** Indeno[1,2,3-cd]pyrene

