

EPA Method 8315A, Carbonyl Compounds

Column: SynergiTM 4 μ m Hydro-RP 80 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4375-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	30	70
	2	20	30	70
	3	35	0	100
	4	50	0	100

Flow Rate: 1.2 mL/min


Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)

Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

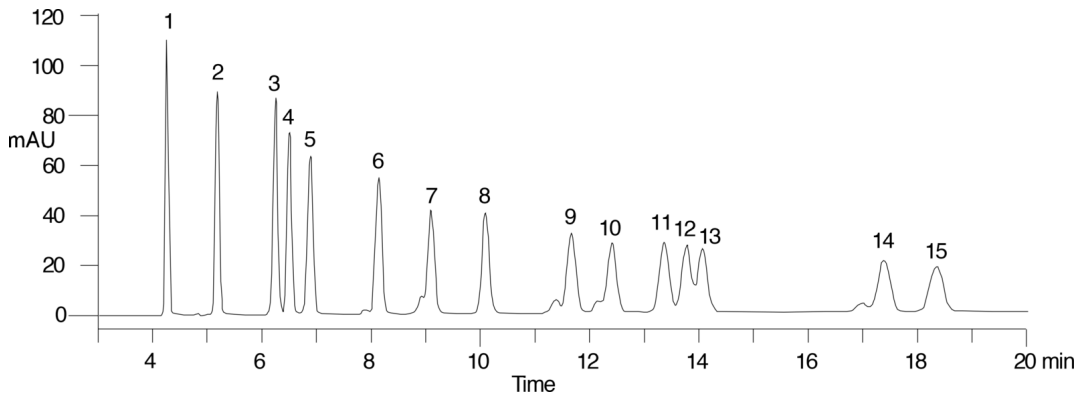

EPA Method 8315A, Carbonyl Compounds

ANALYTES:

- 1** Formaldehyde
- 2** Acetaldehyde
- 3** Acrolein
- 4** Acetone
- 5** Propanal
- 6** Crotonaldehyde
- 7** Butanal
- 8** Benzaldehyde
- 9** Isovaleraldehyde
- 10** Pentanal
- 11** o-Tolualdehyde
- 12** m-Tolualdehyde
- 13** p-Tolualdehyde
- 14** Hexanal
- 15** 2,4-Dimethylbenzaldehyde

