


EPA Method 8015 - Diesel Range Organics

Column: Zebron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G002-11
Oven Profile: 50Å°C for 4 min 300Å°C for 16.5 min 320Å°C for 24.5 min 320Å°C for 30 min
Carrier Gas: Constant Pressure Helium, 1.33 mL/min
Injection: Split 15.0000000000:1 0 µL @ 200°C
Detection: Refractive Index (340°C)

**Products used in this application:****ANALYTES:**

- 1 Diesel fuel
- 2 o-Terphenyl

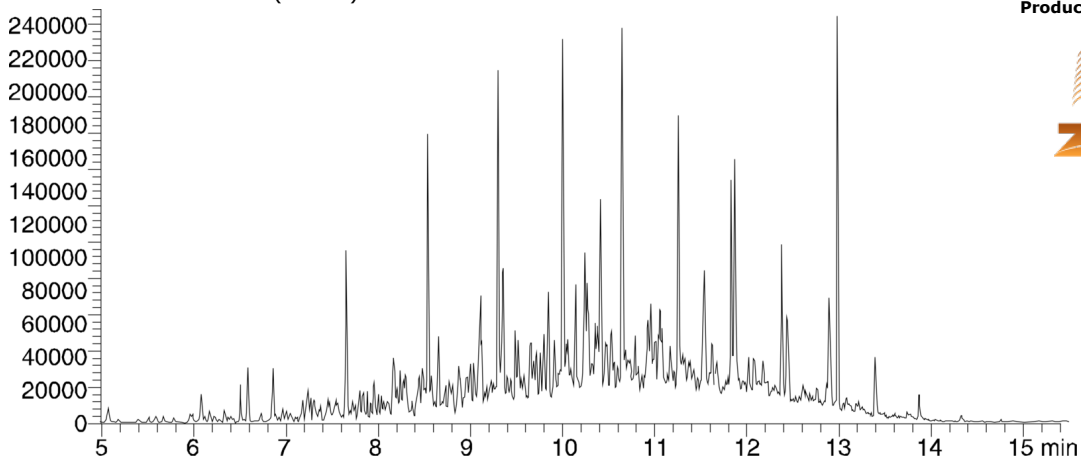