

EPA Method 8330 Explosives mix using Synergi Polar-RP

Column: SynergiTM 4 μ m Polar-RP 80 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4336-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: MeOH

Eluent C: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
	1	0	60	5	35

Flow Rate: 1 mL/min

Col. Temp.: 35 °C

Detection: Refractive Index @ 0.0000000000 (ambient)


Analyst Note: SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, POLAR-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-6076

- Holder Part No.: KJ0-4282



14769

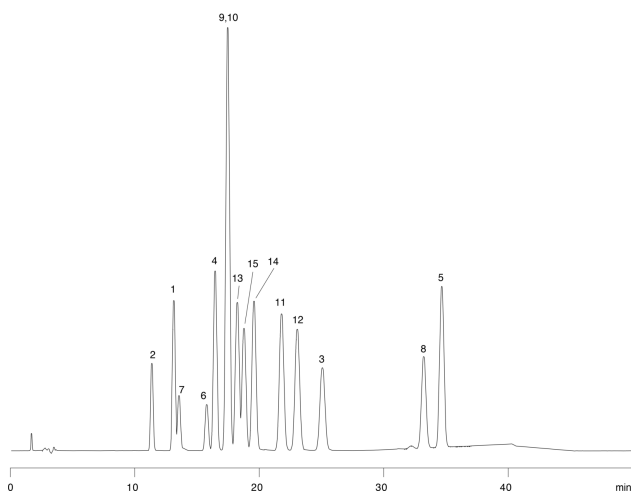

EPA Method 8330 Explosives mix using Synergi Polar-RP

ANALYTES:

- 1** HMX
- 2** RDX
- 3** 1,3,5-Trinitrobenzene
- 4** 1,2-Dinitrobenzene
- 5** Tetryl
- 6** 1,3-Dinitrobenzene
- 7** Nitrobenzene
- 8** 2,4,6-Trinitrotoluene
- 9** 4-Amino-dinitrotoluene
- 10** 2-Amino-dinitrotoluene
- 11** 2,6-Dinitrotoluene
- 12** 2,4-Dinitrotoluene
- 13** 2-Nitrotoluene
- 14** 4-Nitrotoluene
- 15** 3-Nitrotoluene

