


EPA Method 8015: Nonhalogenated Organics on ZB-5

Column: Zebron™ ZB-5, GC Cap. Column 30 m x 0.53 mm x 1.50 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.53 mm x 1.5 µm
Order No: 7HK-G002-28
Oven Profile:

Carrier Gas: Constant Flow Helium, 4.8 mL/min
Injection: Split 20.0000000000:1 0 µL @ 225°C
Detection: Refractive Index (350°C)
14772



Products used in this application:

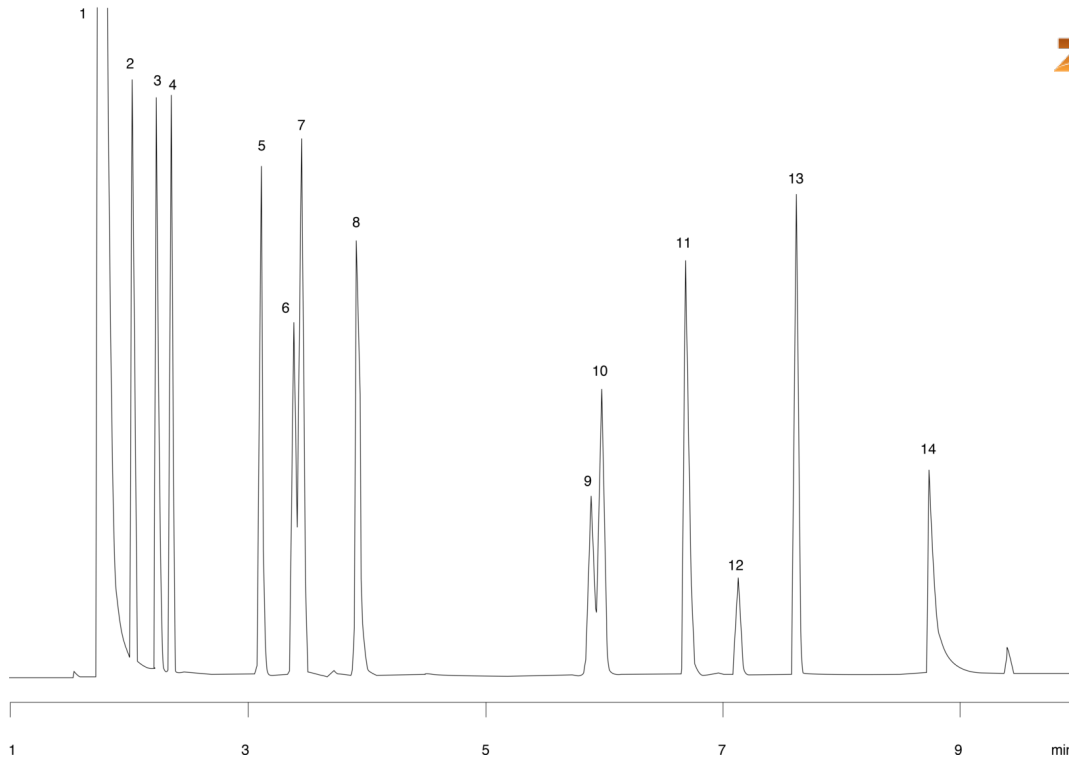

EPA Method 8015: Nonhalogenated Organics on ZB-5

ANALYTES:

- 1 Methanol
- 2 Ethanol
- 3 Acetonitrile
- 4 Diethyl ether
- 5 Propanenitrile
- 6 2-Butanone
- 7 2-Methyl-2-propanitrile
- 8 Isobutanol
- 9 1,4-Dioxane
- 10 Methylmethacrylate
- 11 4-Methyl-2-pentanone
- 12 impurity
- 13 Ethylmethacrylate
- 14 Acrylamide

