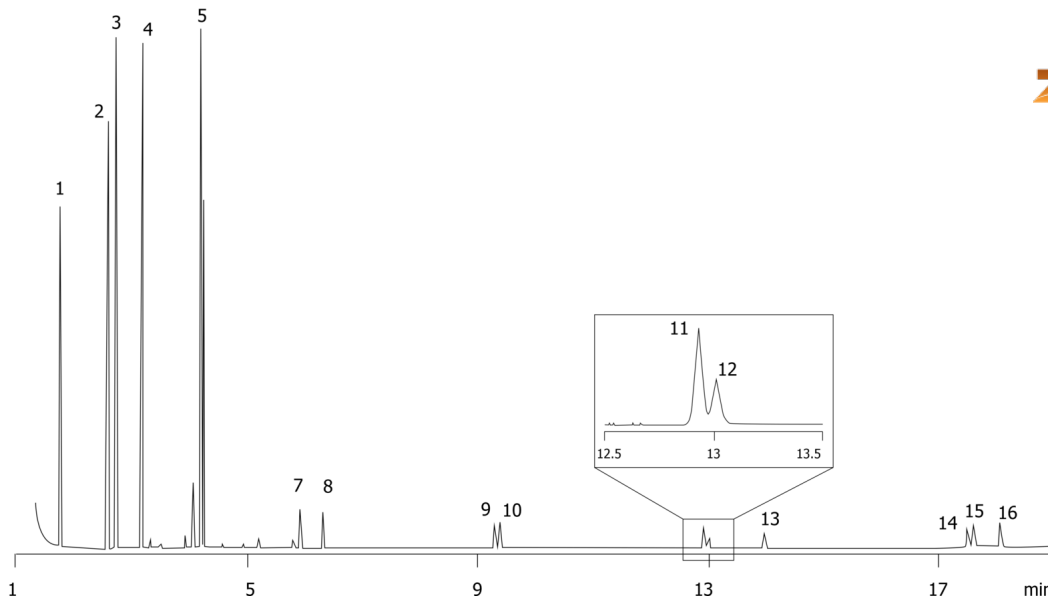


EPA Method 8100: PAH's on ZB-1 30 0.25 0.25**Column:** Zebron™ ZB-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G001-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 20.000000000:1 0 µL @ 260°C**Detection:** Refractive Index (180°C)

14782

**Products used in this application:**

EPA Method 8100: PAH's on ZB-1 30 0.25 0.25

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benz[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Benzo[g,h,i]perylene
- 15 Dibenzo[a,h]anthracene
- 16 Indeno[1,2,3-cd]pyrene

