

PONA: FCC on ZB-1

Column: Zebtron™ ZB-1, GC Cap. Column 100 m x 0.25 mm x 0.50 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 100 meters x 0.25 mm x 0.5 µm

Order No: 7MG-G001-17


Oven Profile:

Carrier Gas: Constant Flow Helium, 1.55 mL/min

Injection: Split 40.0000000000:1 0 µL @ 220°C

Detection: Refractive Index (320°C)

14802

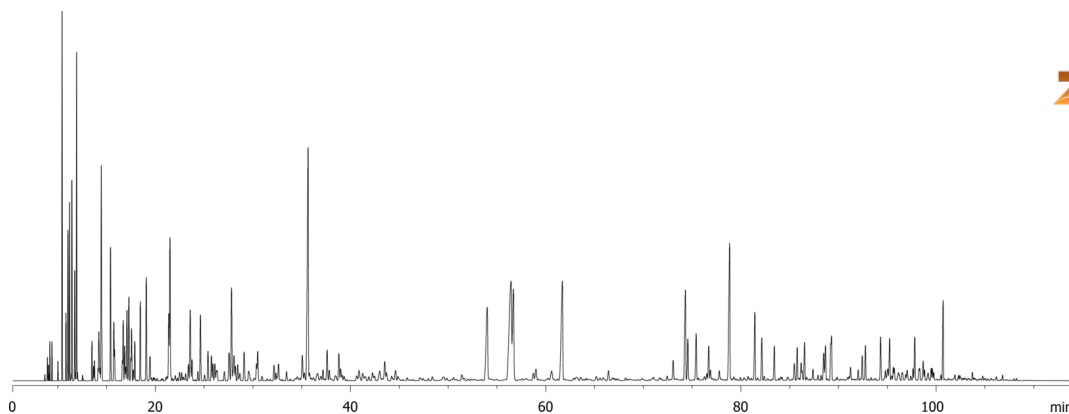

Products used in this application:



PONA: FCC on ZB-1

ANALYTES:

- 1 Pentane
- 2 2,2-Dimethylpentane
- 3 Methylcyclopentane
- 4 Dimethylpentane
- 5 1-Methyl-1-cyclohexene
- 6 Benzene
- 7 Cyclohexane
- 8 2-Methylhexane
- 9 2,3-Dimethylpentane
- 10 3-Methylhexane
- 11 2,2,4-Trimethylpentane
- 12 Heptane
- 13 2,3,4-Trimethylpentane
- 14 Toluene
- 15 2,2,3-Trimethylpentane
- 16 2,3-Dimethylhexane
- 17 Ethylbenzene
- 18 m-Xylene
- 19 p-Xylene
- 20 o-Xylene
- 21 Decane
- 22 Dodecane

