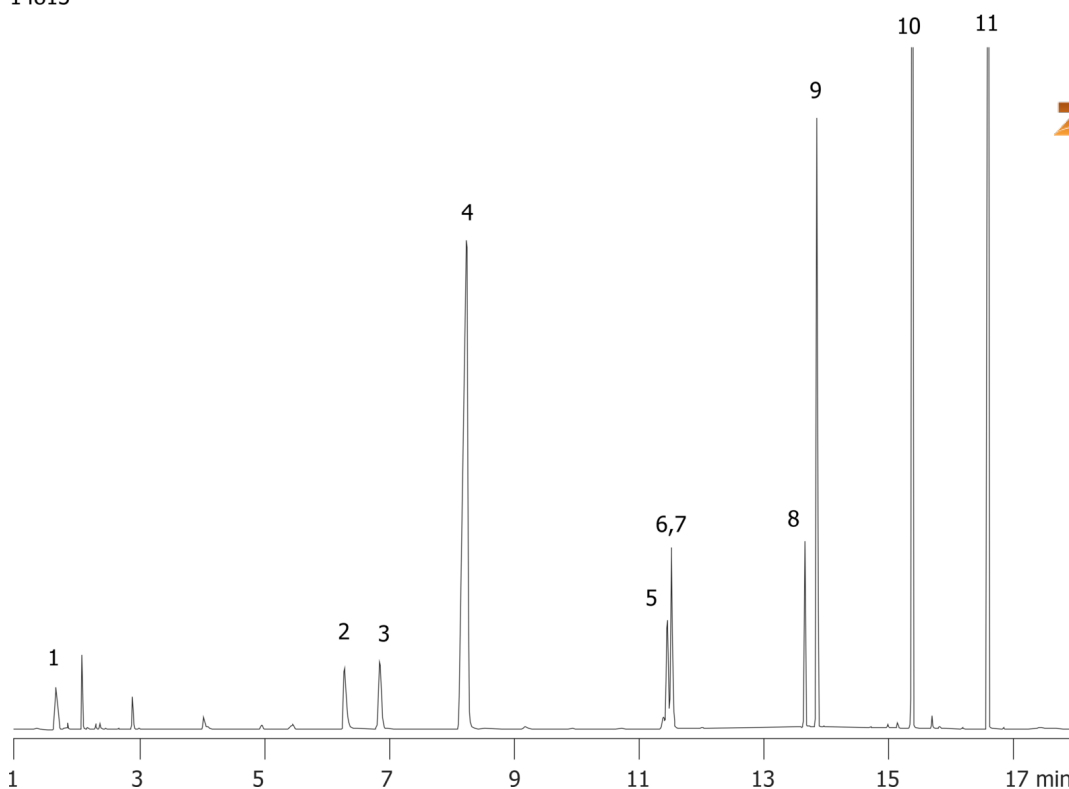


EPA Method 552.2: Haloacetic Acids on ZB-5 30x0.25x0.25**Column:** Zebron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G002-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.9 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 200°C**Detection:** Refractive Index (260°C)

14813

**Products used in this application:**

EPA Method 552.2: Haloacetic Acids on ZB-5 30x0.25x0.25

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 1,2,3-Trichloropropane
- 6 Bromochloroacetic acid
- 7 Trichloroacetic acid
- 8 Dibromoacetic acid
- 9 Bromodichloroacetic acid
- 10 Chlorodibromoacetic acid
- 11 Tribromoacetic acid

