


Nitrosamines - EPA Method 607 on ZB-5 15x0.53x1.5**Column:** Zebron[™] ZB-5, GC Cap. Column 15 m x 0.53 mm x 1.50 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.53 mm x 1.5 µm**Order No:** 7EK-G002-28**Oven Profile:****Carrier Gas:** Constant Flow Helium, 11.4 mL/min**Injection:** Split 15.0000000000:1 0 µL @ 240°C**Detection:** Refractive Index (250°C)

14865

**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosodi-n-propylamine
- 3 N-Nitrosodiphenylamine

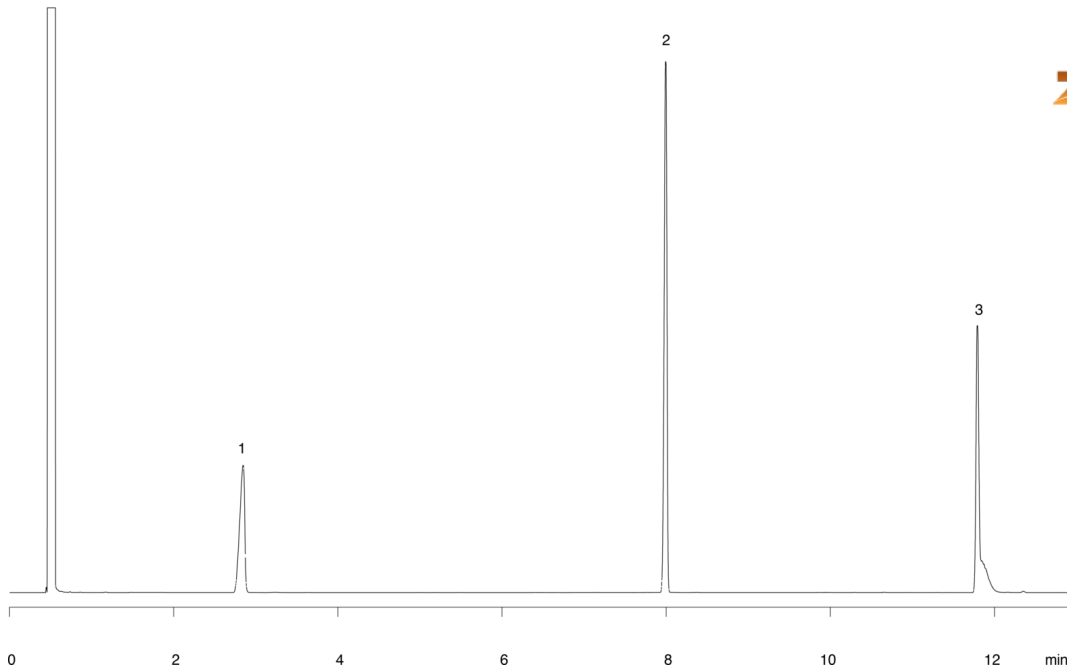