


Nitroaromatics - EPA Method 609: ZB-5 15x0.53x1.5**Column:** ZebronTM ZB-5, GC Cap. Column 15 m x 0.53 mm x 1.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.53 mm x 1.5 μ m**Order No:** 7EK-G002-28**Oven Profile:****Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 15.0000000000:1 0 μ L @ 250°C**Detection:** Refractive Index (300°C)

14866

**Products used in this application:****ANALYTES:**

- 1 Nitrobenzene
- 2 Isophorone
- 3 2,6-Dinitrotoluene
- 4 2,4-Dinitrotoluene

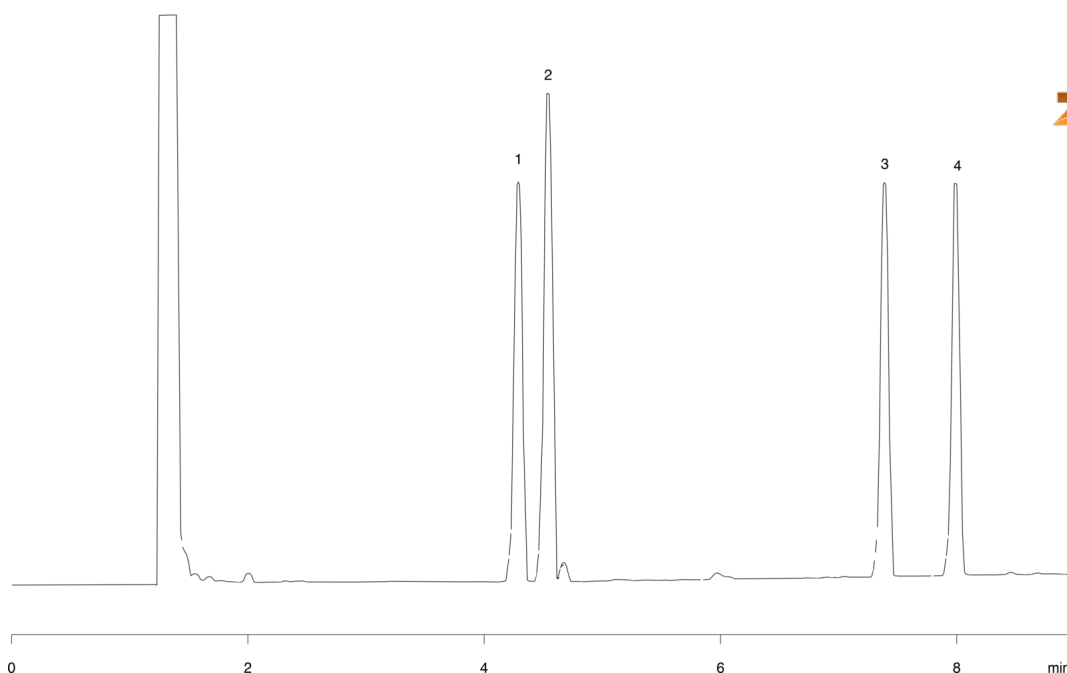