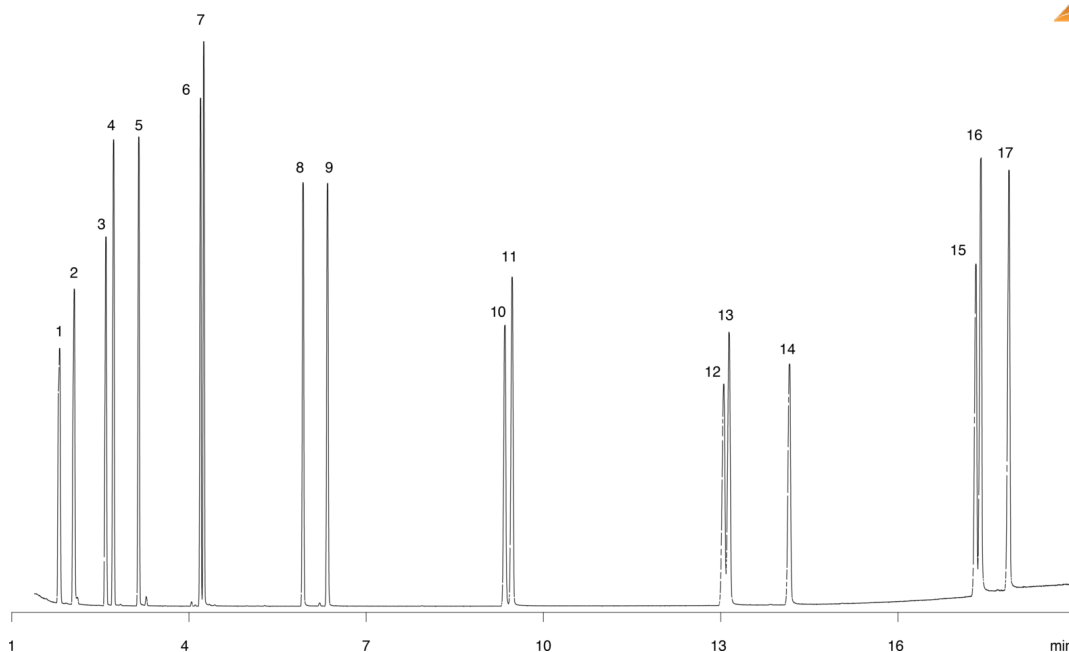


Polyaromatichydrocarbons (PAH):EPA Method 610 - ZB-5 30x0.25x0.25**Column:** Zebtron™ ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G002-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 15.0000000000:1 0 µL @ 260°C**Detection:** Refractive Index (240°C)

14872

**Products used in this application:**

Polyaromatichydrocarbons (PAH):EPA Method 610 - ZB-5 30x0.25x0.25

ANALYTES:

- 1 Naphthalene
- 2 2-Methylnaphthalene
- 3 Acenaphthylene
- 4 Acenaphthene
- 5 Fluorene
- 6 Phenanthrene
- 7 Anthracene
- 8 Fluoranthene
- 9 Pyrene
- 10 Benz[a]anthracene
- 11 Chrysene
- 12 Benzo[b]fluroanthene
- 13 Benzo[k]fluroanthene
- 14 Benzo[a]pyrene
- 15 Indeno[1,2,3-cd]pyrene
- 16 Dibenz[a,h]anthracene
- 17 Benzo[g,h,i]perylene

