


Phthalate Esters: EPA Method 606 - ZB-5 30x0.53x1.5**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.53 mm x 1.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 1.5 μ m**Order No:** 7HK-G002-28**Oven Profile:****Carrier Gas:** Constant Flow Helium, 12.9 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 300°C**Detection:** Refractive Index (300°C)

14876

**Products used in this application:****ANALYTES:**

- 1 Dimethyl phthalate
- 2 Diethyl phthalate
- 3 Di-n-butyl phthalate
- 4 Butyl benzyl phthalate
- 5 bis(2-Ethylhexyl)phthalate
- 6 Di-n-octyl phthalate

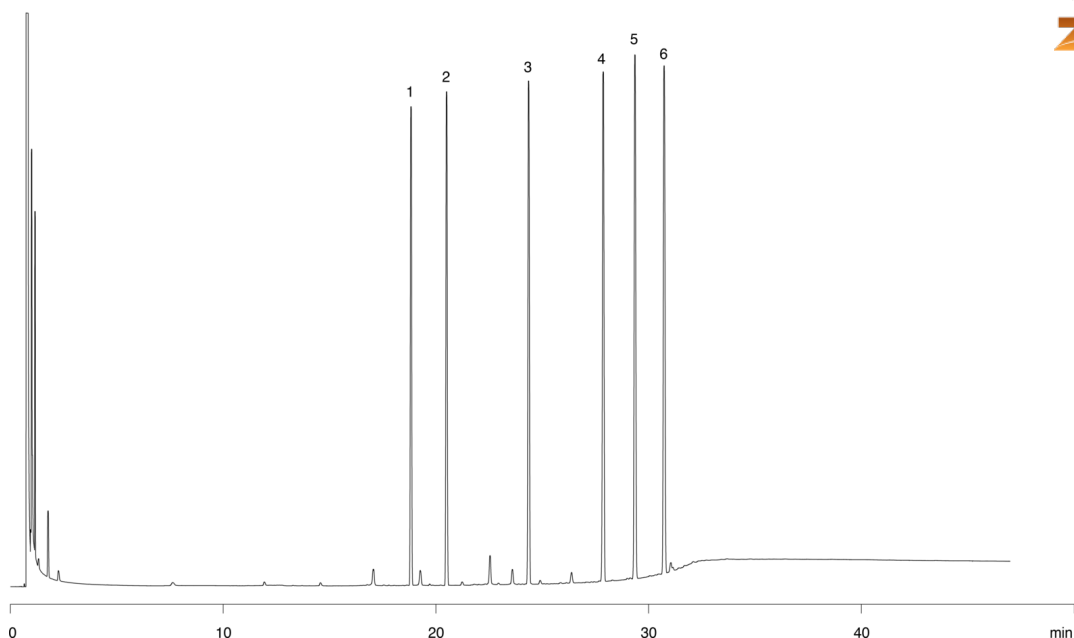