

## Triazine Herbicides: EPA Method 619 - ZB-5 30x0.32x0.5

**Column:** Zebtron™ ZB-5, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea

**Phase:** 5% Phenyl 95% Dimethylpolysiloxane

**Dimensions:** 30 meters x 0.32 mm x 0.5 µm

**Order No:** 7HM-G002-17

**Oven Profile:**

**Carrier Gas:** Constant Flow Helium, 1.2 mL/min


**Injection:** Split 20.0000000000:1 0 µL @ 250°C

**Detection:** Refractive Index (240°C)

14885



Products used in this application:

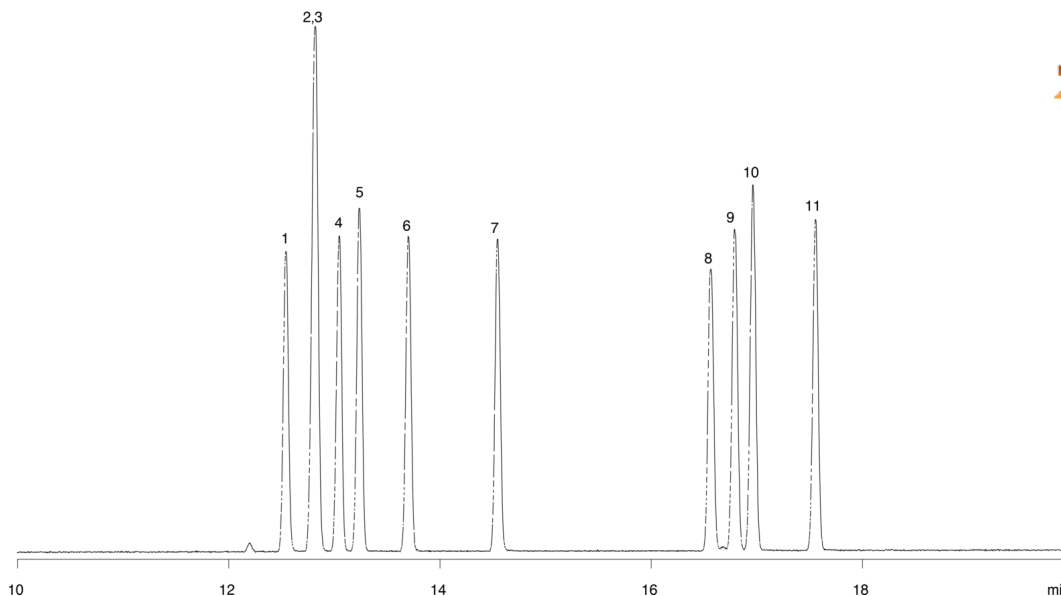

## Triazine Herbicides: EPA Method 619 - ZB-5 30x0.32x0.5

### ANALYTES:

- 1 Atraton
- 2 Simazine
- 3 Prometon
- 4 Atrazine
- 5 Propazine
- 6 Terbutylazine
- 7 Secbumetone
- 8 Simetryn
- 9 Ametryn
- 10 Prometryn
- 11 Terbutryn

