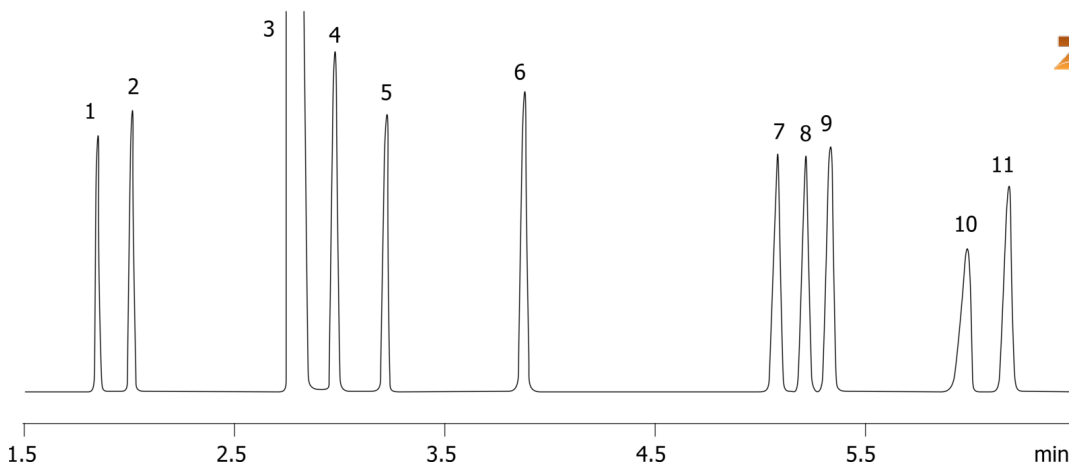


BTEX in Under Six Minutes on ZB-WAX**Column:** ZebtronTM ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G007-17**Oven Profile:** 75 $^{\circ}$ C for 1 min 90 $^{\circ}$ C for 6 min 90 $^{\circ}$ C for 9 min**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 250 $^{\circ}$ C**Detection:** Refractive Index (300 $^{\circ}$ C)

14895

**Products used in this application:**

BTEX in Under Six Minutes on ZB-WAX

ANALYTES:

- 1 Pentane
- 2 Heptane
- 3 Methylene chloride
- 4 Benzene
- 5 Decane
- 6 Toluene
- 7 Ethylbenzene
- 8 p-Xylene
- 9 m-Xylene
- 10 Dodecane
- 11 o-Xylene

