


Phenols: EPA Method 604 - ZB-1701 30x0.25x0.25**Column:** ZebronTM ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G006-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 300°C**Detection:** Refractive Index (300°C)

14897

**Products used in this application:****ANALYTES:**

- 1 2-Chlorophenol
- 2 Phenol
- 3 2-Nitrophenol
- 4 2,4-Dimethylphenol
- 5 2,4-Dichlorophenol
- 6 2,4,6-Trichlorophenol
- 7 4-Chloro-3-methylphenol
- 8 2-Methyl-4,6-dinitrophenol
- 9 4-Nitrophenol

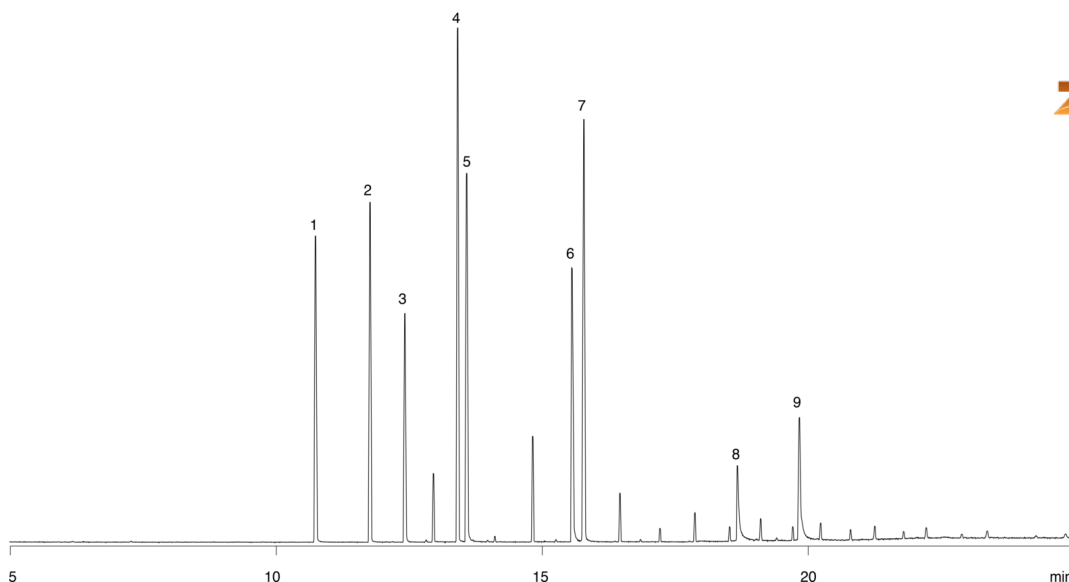