


Cocaine and Barbiturates on ZB-5**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G002-11**Oven Profile:** 75°C to 320°C at 25°C/min.**Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Splitless :1 2 µL @ 250°C**Detection:** Mass Selective (MSD) (180°C)

14904

**Products used in this application:****ANALYTES:**

- 1 Pentobarbital
- 2 Meprobamate
- 3 Hexobarbital
- 4 Cocaine
- 5 Diazepam
- 6 6-MAM

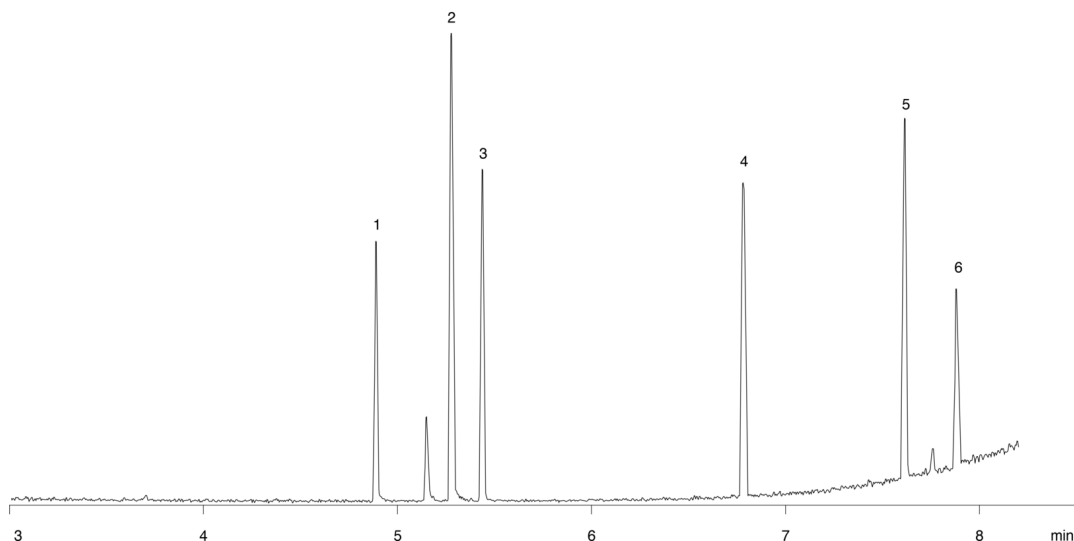