

## Jupiter C18 P4\_2

**Column:** Jupiter® 5 µm C18 300 Å, LC Column 150 x 2 mm, Ea

**Dimensions:** 150 x 2 mm ID

**Order No:** 00F-4053-B0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA in 95:5 Water/ACN

**Eluent B:** 0.085% TFA in 75:20:5 ACN/IPA/Water

| Gradient        | Step No. | Time (min) | Pct A | Pct B |
|-----------------|----------|------------|-------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 80    | 20    |
|                 | <b>2</b> | 20         | 20    | 80    |

**Flow Rate:** 0.22 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C18 4 x 2.0mm, 10/Pk Part No.: AJ0-4320


- Holder Part No.: KJ0-4282



Products used in this application:



14929

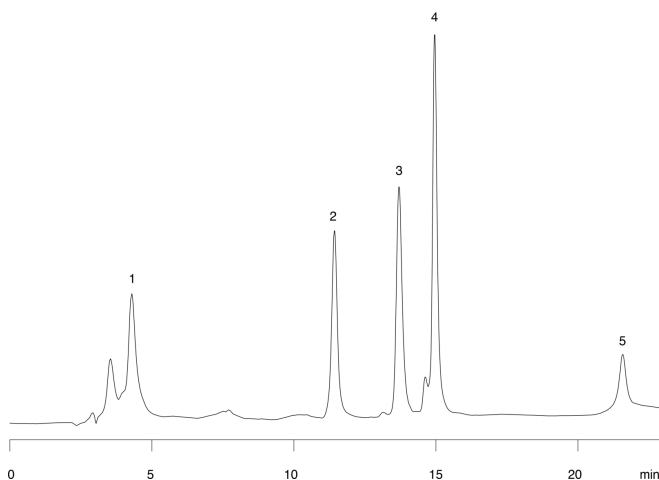

### ANALYTES:

- 1** Aprotinin  
Retention Time: 4.29 min
- 2** Ribonuclease  
Retention Time: 11.43 min
- 3** Apotransferrin  
Retention Time: 13.7 min
- 4** Lactalbumin  
Retention Time: 14.95 min
- 5** Actin  
Retention Time: 21.6 min

