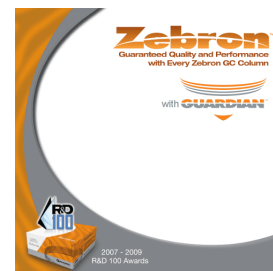
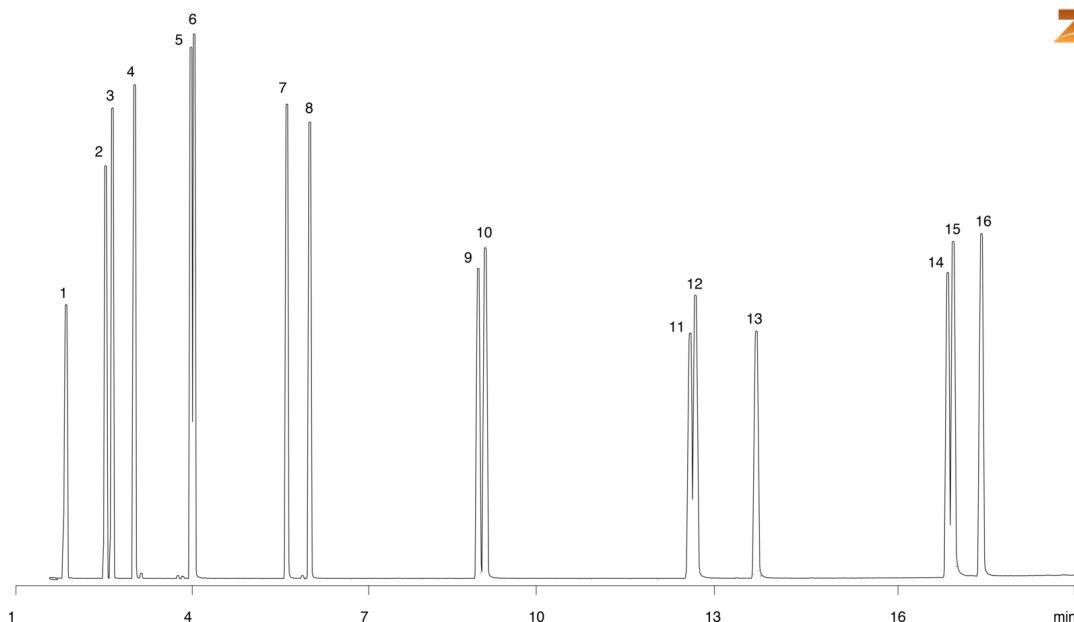


Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-5**Column:** Zebtron[™] ZB-5 w/GUARDIAN[™], 5 m, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G002-11-GGA**Oven Profile:** 240 $^{\circ}$ C for 4.67 min 275 $^{\circ}$ C for 13.42 min 320 $^{\circ}$ C for 17.92 min 320 $^{\circ}$ C for 22.92 min**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 15.0000000000:1 0 μ L @ 260 $^{\circ}$ C**Detection:** Refractive Index (280 $^{\circ}$ C)

14942

**Products used in this application:**

Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-5

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benz[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Indeno[1,2,3-cd]pyrene
- 15 Dibenzo[a,h]anthracene
- 16 Benzo[g,h,i]perylene

