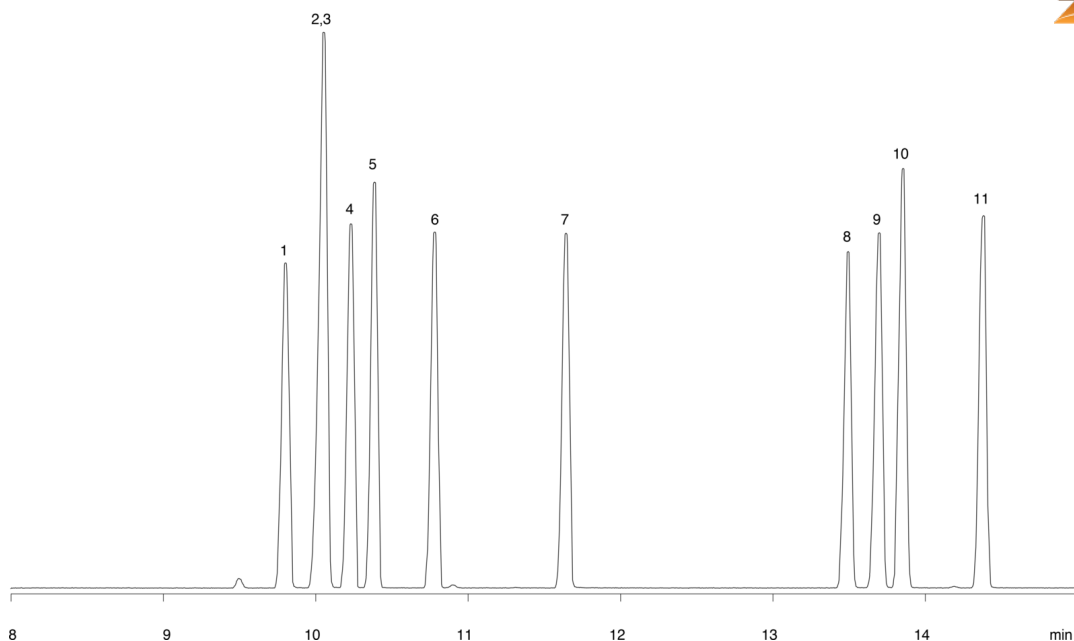


Triazine Herbicides on ZB-5MS**Column:** Zebtron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:** 220Å°C for 20 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (240°C)

14954

**Products used in this application:**

Triazine Herbicides on ZB-5MS

ANALYTES:

- 1 Atraton
- 2 Simazine
- 3 Prometon
- 4 Atrazine
- 5 Propazine
- 6 Terbutylazine
- 7 Secbumetone
- 8 Simetryn
- 9 Ametryn
- 10 Prometryn
- 11 Terbutryn

