


Spearmint Oil on ZB-5MS

Column: Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G010-11
Oven Profile: 40Å°C for 4 min 120Å°C for 14 min 200Å°C for 19.33 min 300Å°C for 23.33 min
Carrier Gas: Constant Flow Helium, 0.8 mL/min
Injection: Split 25.0000000000:1 0 µL @ 250°C
Detection: Refractive Index (240°C)

14958



Products used in this application:

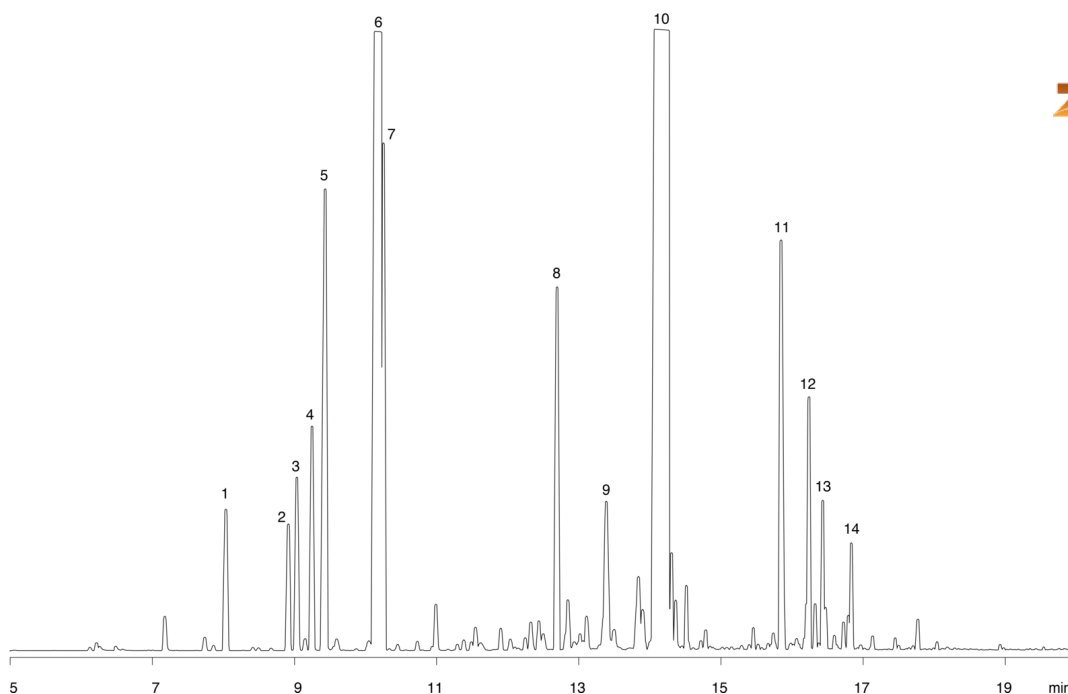

Spearmint Oil on ZB-5MS

ANALYTES:

- 1 α -Pinene
- 2 Sabinene
- 3 beta-Pinene
- 4 beta-Myrcene
- 5 3-Octanol
- 6 Limonene
- 7 Eucalyptol
- 8 Menthone
- 9 Dihydrocarvone
- 10 Carvone
- 11 Bourbonene
- 12 Caryophyllene
- 13 Farnesene
- 14 Germacrene

