

Estradiol - USP Method

Column: Gemini® 5 µm C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	45	55

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: Refractive Index @ 0.0000000000 (25 °C)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7597

- Holder Part No.: KJ0-4282



Products used in this application:

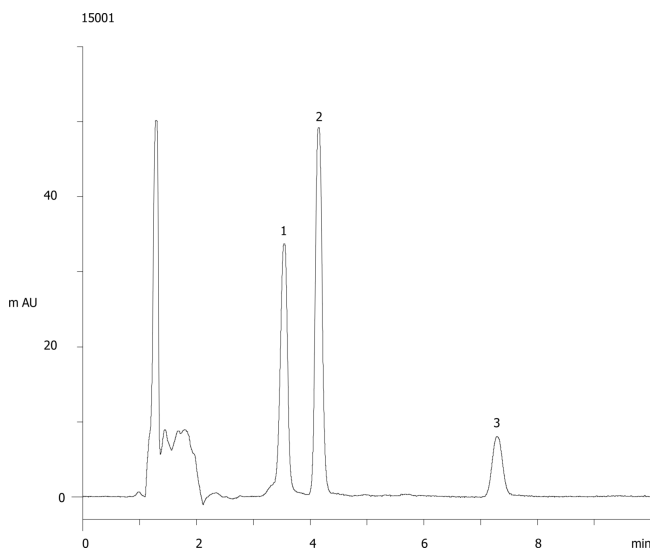

ANALYTES:

1 Ethyl paraben

Retention Time: 3.542 min

2 Estrone

Retention Time: 4.149 min

3 Estradiol

Retention Time: 7.293 min

