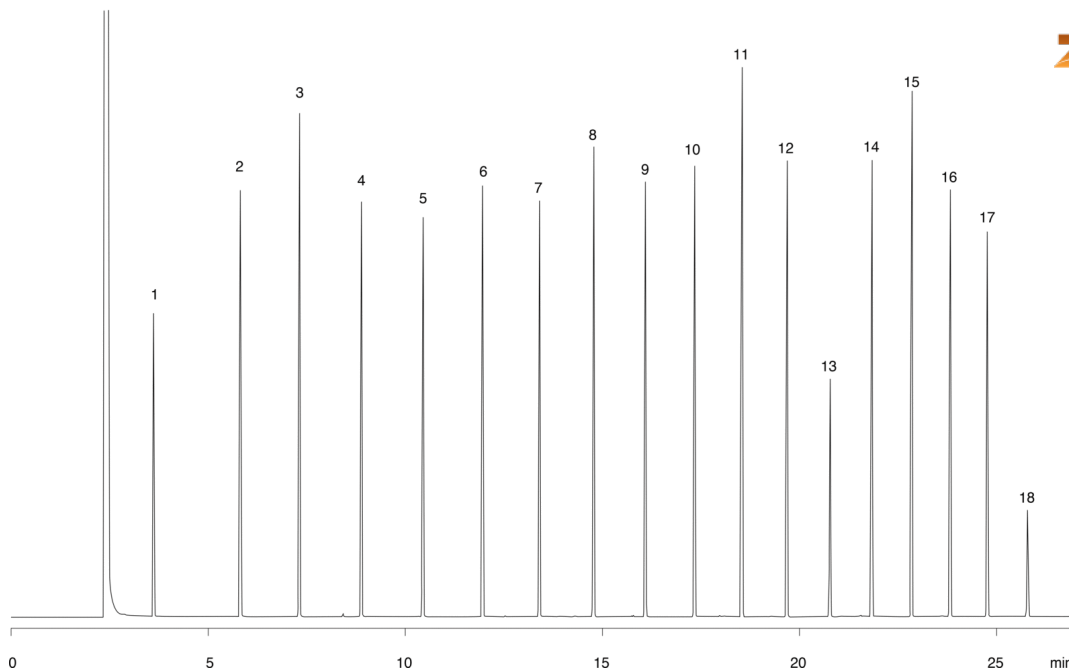


FAMES on ZB-5MS**Column:** ZebtronTM ZB-5MS, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 280°C**Detection:** Refractive Index (325°C)

15423

**Products used in this application:**

FAMES on ZB-5MS

ANALYTES:

- 1 Methyl caproate (C6)
- 2 Methyl caprylate (C8)
- 3 Methyl nonoate (C9)
- 4 Methyl caprate (C10)
- 5 Methyl undecanoate (C11)
- 6 Methyl laurate (C12)
- 7 Methyl tridecanoate (C13)
- 8 Methyl myristate (C14)
- 9 Methyl pentadecanoate (C15)
- 10 Methyl palmitate (C16)
- 11 Methyl heptadecanoate (C17)
- 12 Methyl stearate (C18)
- 13 Methyl nonadecanoate (C19)
- 14 Methyl arachidate (C20)
- 15 Methyl heneicosanoate (C21)
- 16 Methyl behenate (C22)
- 17 Methyl tricosanoate (C23)
- 18 Methyl lignocerate (C24)

