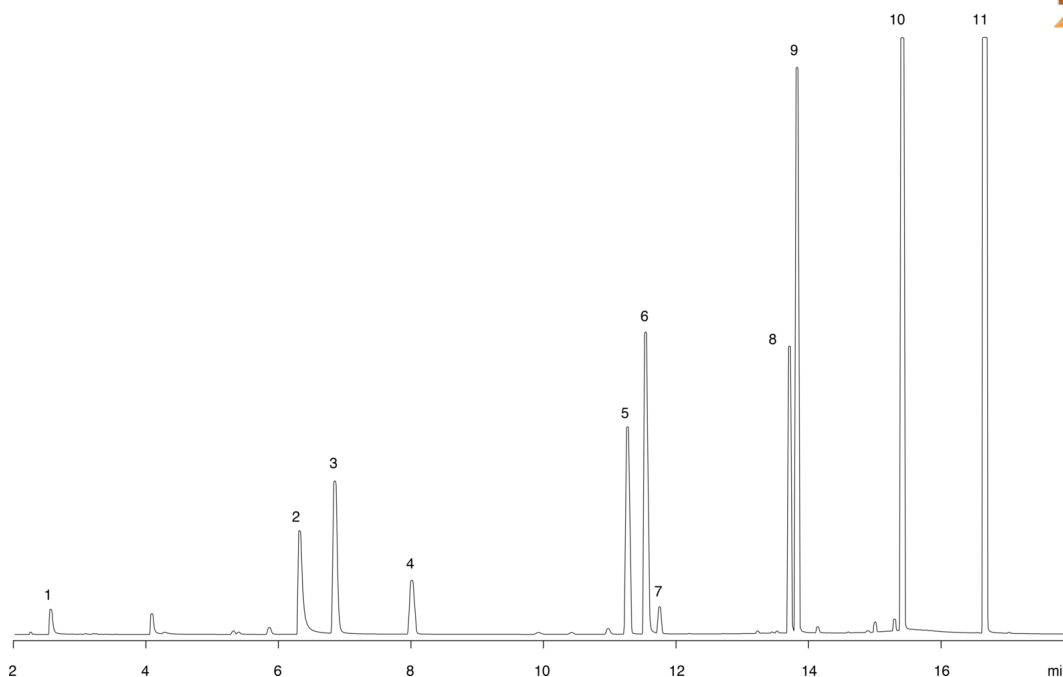


Haloacetic Acids - EPA Method 552.2**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G010-11**Oven Profile:** 45 $^{\circ}$ C for 10 min 145 $^{\circ}$ C for 16.67 min 145 $^{\circ}$ C for 21.67 min**Carrier Gas:** Constant Flow Helium, 1.9 mL/min**Injection:** Split 20.0000000000:1 0 μ L @ 200 $^{\circ}$ C**Detection:** Refractive Index (260 $^{\circ}$ C)

15431

**Products used in this application:**

Haloacetic Acids - EPA Method 552.2

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 Bromochloroacetic acid
- 7 1,2,3-Trichloropropane
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 Tribromoacetic acid

