


Aliphatic Amines on ZB-5MS

Column: Zebron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 15.0000000000:1 0 µL @ 220°C**Detection:** Refractive Index (240°C)

15432



Products used in this application:

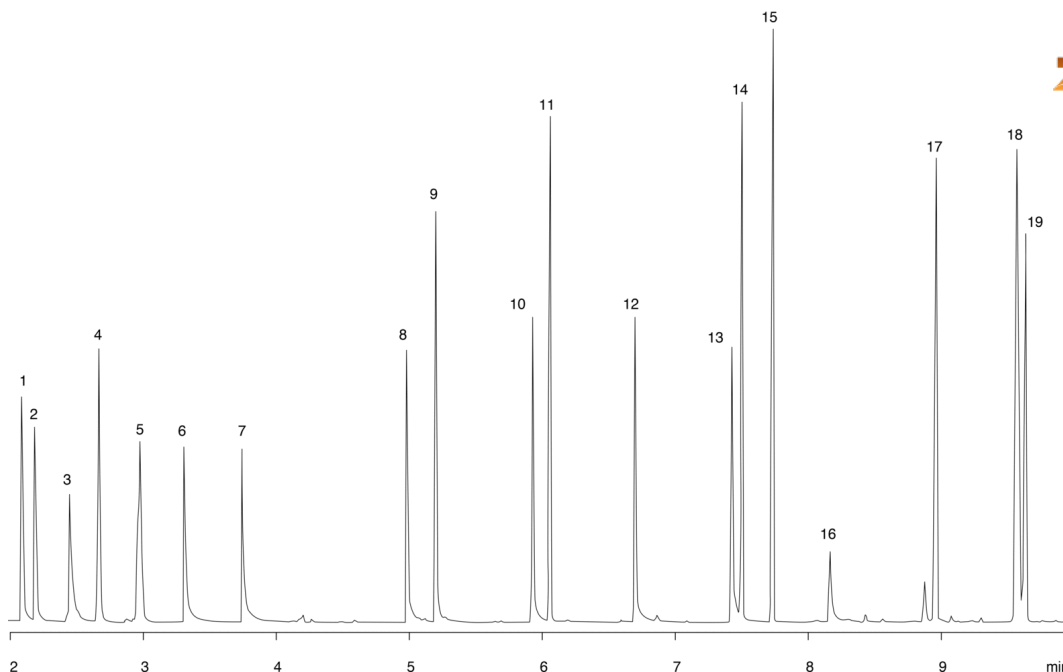

Aliphatic Amines on ZB-5MS

ANALYTES:

- 1 sec-Butylamine
- 2 Isobutylamine
- 3 n-Butylamine
- 4 Diisopropylamine
- 5 Triethylamine
- 6 Isopentylamine
- 7 Pentylamine
- 8 Hexylamine
- 9 Cyclohexylamine
- 10 Heptylamine
- 11 Dibutylamine
- 12 Octylamine
- 13 Nonylamine
- 14 Dipentylamine
- 15 Tributylamine
- 16 Decylamine
- 17 Dihexylamine
- 18 Dicyclohexylamine
- 19 contaminant

