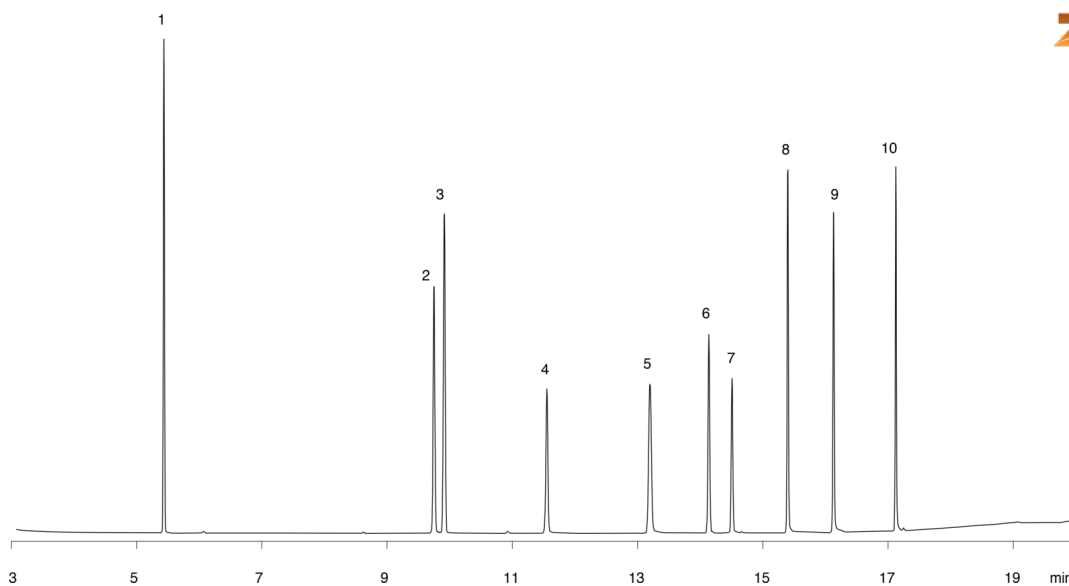


Nitrogen and Phosphorous Pesticides on ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (240°C)

15440

**Products used in this application:**

Nitrogen and Phosphorous Pesticides on ZB-5MS

ANALYTES:

- 1 Vernolate
- 2 Terbufos
- 3 Diazinon
- 4 Metribuzin
- 5 Metolachlor
- 6 MGK-264
- 7 MGK-264
- 8 Butachlor
- 9 Carboxin
- 10 Norflurazon

