


Nitrogen and Phosphorous Pesticides on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (240°C)

15442

**Products used in this application:****ANALYTES:**

- 1 Propachlor
- 2 Trifluralin
- 3 Benefin
- 4 Profluralin
- 5 Isopropalin
- 6 Pendimethalin
- 7 Oxadiazon
- 8 Oxyfluorfen

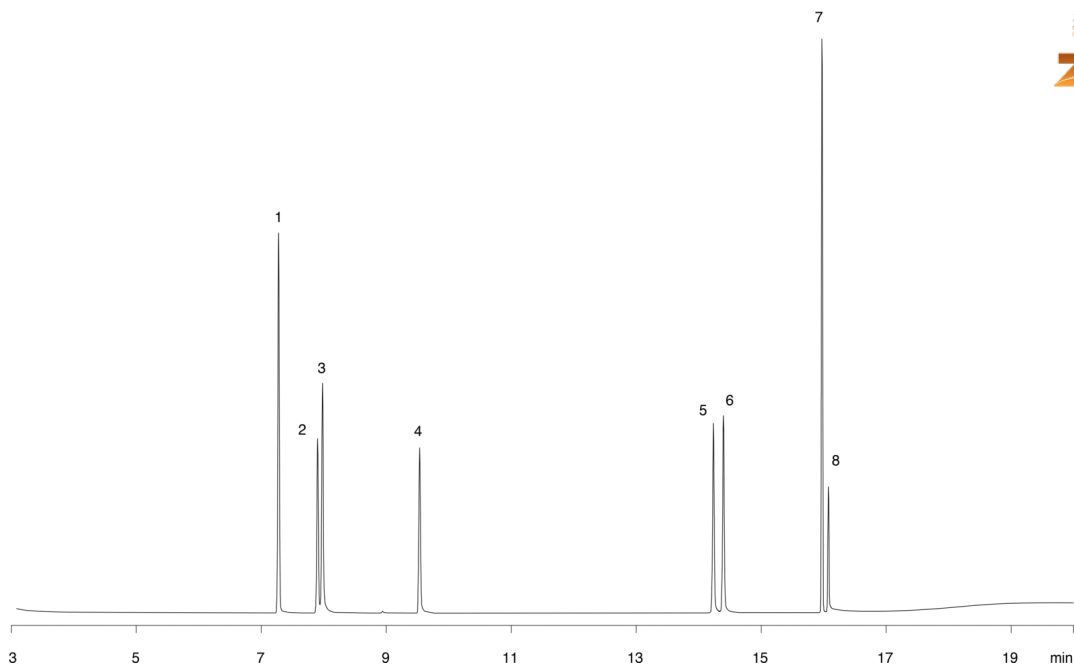